

Europäisches Patentamt **European Patent Office** Office européen des brevets



EP 1 442 715 A2 (11)

(12)

# **EUROPEAN PATENT APPLICATION**

(43) Date of publication: 04.08.2004 Bulletin 2004/32

(21) Application number: 04001428.4

(22) Date of filing: 23.01.2004

(51) Int Cl.7: **A61B 17/17**, A61B 19/00, A61F 2/32, A61F 2/44, A61F 2/46

(84) Designated Contracting States: AT BE BG CH CY CZ DE DK EE ES FI FR GB GR HU IE IT LI LU MC NL PT RO SE SI SK TR **Designated Extension States:** AL LT LV MK

(30) Priority: **30.01.2003 US 354562** 

25.04.2003 US 423515

(71) Applicant: SURGICAL NAVIGATION TECHNOLOGIES, INC. Louisville, Colorado 80027 (US)


(72) Inventor: Hunter, Mark W. **Broomfield Colorado 80020 (US)** 

(74) Representative: Edlund, Fabian Awapatent AB, P.O. Box 11 394 404 28 Göteborg (SE)

#### Tunable spinal implant and apparatus for its post-operative tuning (54)

(57)A tunable implant, system, and method enables a tunable implant (302) to be adjusted within a patient. The tunable implant includes a securing mechanism (308,310) to secure the implant in the patient, a actuation portion that enables the implant to move and an adjustment portion that permits adjustment of the implant after the implant has been positioned within the patient. The method of adjusting the tunable implant includes analyzing the operation of the implant, determining if

any adjustments are necessary and adjusting the implant to improve implant performance. The implant system includes both the tunable implant and a telemetric system (378,384) that is operable to telemetrically receive data from the tunable implant where the data is used to determine if adjustment of the tunable implant is necessary. The system also includes an instrument assembly that is used for performing spinal surgery where the instrument assembly includes a mounting platform (268) and a jig (270).

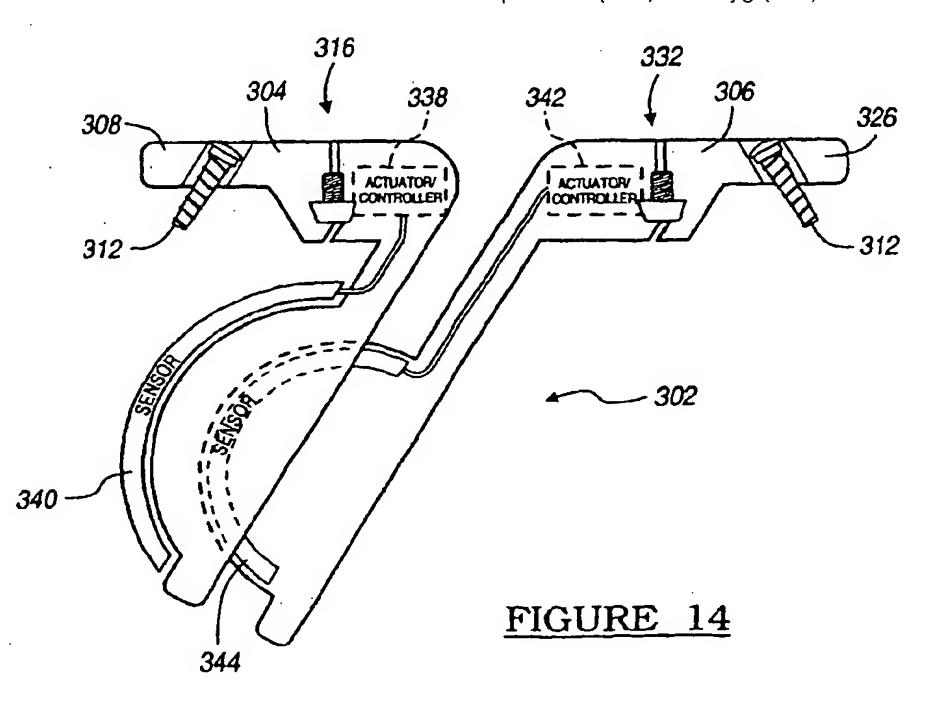

Printed by Jouve, 75001 PARIS (FR)

# **Description**

### FIELD OF THE INVENTION

[0001] The present invention generally relates to adjusting or tuning an implant, and more specifically, to pre-operative planning to select an implant and technique and post-operative tuning of the implant by a telemetric or a minimally invasive procedure.

#### BACKGROUND OF THE INVENTION

Image guided medical and surgical procedures utilize patient images obtained prior to or during a medical procedure to guide a physician performing the procedure. Recent advances in imaging technology, especially in imaging technologies that produce highly-detailed, computer-generated two, three and four-dimensional images, such as computed tomography (CT), magnetic resonance imaging (MRI), isocentric C-arm fluoroscopic imaging, fluoroscopes or ultrasounds have increased the interest in image guided medical procedures. During these image guided medical procedures, the area of interest of the patient that has been imaged is displayed on a display. Surgical instruments and/or implants that are used during this medical procedure are tracked and superimposed onto this display to show the location of the surgical instrument relative to the area of interest in the body.

[0003] Other types of navigation systems operate as an image-less system, where an image of the body is not captured by an imaging device prior to the medical procedure. With this type of procedure, the system may use a probe to contact certain landmarks in the body, such as landmarks on bone, where the system generates either a two-dimensional or three-dimensional model of the area of interest based upon these contacts. This way, when the surgical instrument or other object is tracked relative to this area, they can be superimposed on this model.

[0004] Most types of orthopedic medical procedures are performed using conventional surgical techniques. These techniques generally involve opening the patient in a relatively invasive manner to provide adequate viewing by the surgeon during the medical procedure. These types of procedures, however, generally extend the recovery period for the patient due to the extent of soft tissue and muscular incisions resulting from the medical procedure. Use of image guided technology in orthopedic medical procedures would enable a more minimally invasive type of procedure to be performed to thereby reduce the overall recovery time and cost of the procedure. Use of the image guided procedure may also enable more precise and accurate placement of an implant within the patient.

[0005] Once the implant has been surgically positioned within the patient, the patient's surrounding anatomy generally heals over time with the surrounding skel-

etal and muscular structure regaining a healthy state. However, since the implant is generally implanted when the patient is dysfunctional, this muscular and skeletal adjustment or healing may effect the subsequent range of motion, effectiveness, life expectancy of the implant, performance of the implant, and potentially cause deterioration of surround discs or implants. For example, in a spinal implant, upon the abdominal and back muscles strengthening after the implant procedure, the spine may subsequently align. This alignment may result in the implant or articulation faces of the implant being impinged because of the resultant alignment. This may result in a revision-type surgery that requires the implant to be removed and a subsequent implant being repositioned at the implant site.

[0006] The surgical procedures performed during orthopedic medical procedures, including spinal procedures, require the use of various instruments, assemblies and jigs to perform the procedure. Typically, jigs are used to support a single instrument that must be attached to the area of interest when the instrument is being used. Multiple jigs are thus typically required to be attached and removed from the area of interest as the procedure progresses. Use of multiple jigs and instruments, along with attaching and reattaching to the area of interest provides for a tedious and time consuming procedure. Moreover, inherent inaccuracies due to this procedure may provide less than acceptable results.

[0007] It is, therefore, desirable to provide a method and apparatus for post-operative adjustment or tuning of an implant, such as a spinal implant using telemetric or minimally invasive techniques. It is also desirable to provide an instrument assembly that may be attached to the implant site, such as a spinal implant site, once during the entire procedure, thereby reducing surgical time, costs, as well as increasing surgical accuracy. It is, therefore, an object of the present invention to provide such methods and apparatus for use in medical procedures.

# SUMMARY OF THE INVENTION

[0008] In accordance with the teachings of the present invention, a tunable implant, system, method and associated instruments for use in implanting and adjusting the tunable implant after the implant has been positioned within the patient is disclosed. The tunable implant may be any type of implant, such as a spinal implant.

[0009] In one embodiment, a method for tuning an implant positioned within a patient is provided. This method includes analyzing the operation of the implant that is positioned within the patient, determining if any adjustment of the implant is necessary, and adjusting the implant with the implant positioned within the patient to improve implant performance.

[0010] In another embodiment, a tunable implant is positioned within the patient. The tunable implant in-

Sist &

40

cludes a securing mechanism that is used to secure the implant within the patient. An actuation portion is used to permit the tunable implant to move. An adjustment portion is used to permit adjustment of the tunable implant after the tunable implant is positioned within the patient.

[0011] In yet another embodiment, a tunable implant system for use in adjusting a tunable implant of a patient is provided. This tunable implant system includes a tunable implant having an adjustable portion that is operable to permit adjustment of the tunable implant after the implant is positioned within the patient. A telemetric system is provided and operable to telemetrically receive data from the tunable implant where the data is used to determine adjustment of the tunable implant.

[0012] Still another embodiment, an instrument assembly for use in performing spinal surgery is provided. This instrument includes a mounting platform operable to be positioned adjacent to vertebrae. The instrument assembly also includes a jig that is operable to be removably attached to the mounting platform. The jig is also operable to support an instrument used during the spinal surgery.

[0013] The invention also relates to a display, a navigation system and a method for guiding a medical device to a target in a patient during a medical procedure, also new to the art.

[0014] Further areas of applicability of the present invention will become apparent from the detailed description provided hereinafter. It should be understood that the detailed description and specific examples, while indicating the preferred embodiment of the invention, are intended for purposes of illustration only and are not intended to limit the scope of the invention.

## BRIEF DESCRIPTION OF THE DRAWINGS

[0015] The present invention will become more fully understood from the detailed description and the accompanying drawings, wherein:

[0016] Figure 1 is a diagram of a navigation system employing a display according to the teachings of the present invention;

[0017] Figures 2a and 2b are diagrams representing fundistorted and distorted views of a fluoroscopic C-arm imaging device;

[0018] Figures 3a and 3b is a logic block diagram illustrating a method for employing the display according to the teachings of the present invention;

[0019] Figures 4a-4e illustrate a medical procedure employing the display according to the teachings of the present invention;

[0020] Figure 5 is a figure of the display according to the teachings of the present invention;

[0021] Figure 6 is a split screen view of the display according to the teachings of the present invention;

[0022] Figure 7 is an additional split screen view of the display according to the teachings of the present in-

vention;

[0023] Figures 8a-8g illustrate another medical procedure employing the display according to the teachings of the present invention;

[0024] Figure 9 is an illustration of a dual display according to the teachings of the present invention;

[0025] Figure 10 is a logic block diagram illustrating a method for pre-operative planning and post-operative exam and tuning of an implant according to the teachings of the present invention:

[0026] Figure 11 is a perspective view of a platform and jig used in a minimally invasive surgical navigation spinal procedure;

[0027] Figure 12 is a perspective view of the jig that is operable to be attached to the platform of Figure 11; [0028] Figure 13 is a side view of a cervical disc implant having a minimally invasive adjustment mechanism;

[0029] Figure 14 is a side view of a cervical disc implant having a telemetric adjustment mechanism;

[0030] Figure 15 illustrates an implant of Figures 14 and 15 implanted into a spine;

[0031] Figure 16 is a side cross-sectional view of a cervical disc implant according to the teachings of the present invention;

[0032] Figures 16a-16c are an unfolded and partially folded view of other embodiment of a cervical disc implant according to the teachings of the present invention;

[0033] Figure 17 is a side cross-sectional view of another cervical disc implant according to the teachings of the present invention;

[0034] Figure 18 illustrates a cervical disc system employing multiple cervical disc implants according to the teachings of the present invention;

[0035] Figure 19 illustrates a transmit/receive module used during the motion analysis study of a patient according to the teachings of the present invention; and [0036] Figure 20 illustrates a home based transmit/receive module used for a motion analysis study according to the teachings of the present invention.

# DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENTS

[0037] The following description of the embodiment (s) is merely exemplary in nature and is in no way intended to limit the invention, its application, or uses. Moreover, while the invention is discussed in detail below in regard to orthopedic/spinal surgical procedures, the present invention may be used with any type of medical procedure, including orthopedic, cardiovascular, neurovascular, soft tissue procedures, or any other medical procedures.

[0038] Figure 1 is a diagram illustrating a five or six degree of freedom (5 or 6 DOF) alignment display 10 employed with an image guided navigation system 12 for use in navigating a surgical instrument or implant

during a medical procedure. It should also be noted that the display 10 may be used or employed in an imageless based navigation system, further discussed herein. The navigation system 12 may be used to navigate any type of instrument or delivery system, such as a reamer, impactor, cutting block, saw blade, catheter, guide wires, needles, drug delivery systems, and cell delivery systems. The navigation system 12 may also be used to navigate any type of implant including orthopedic implants, spinal implants, cardiovascular implants, neurovascular implants, soft tissue implants, or any other devices implanted in a patient 14. The navigation system 12 may also be used to navigate implants or devices that are formed as an assembly or from multiple components where the location and orientation of each component is dependent upon one another to be effective in its use. For example, during a spinal procedure, the display may be used to track and align a spinal screw with a spinal rod to insure attachment of each device.

[0039] The navigation system 12 includes an imaging device 16 that is used to acquire pre-operative or realtime images of the patient 14. The imaging device 16 is a fluoroscopic C-arm x-ray imaging device that includes a C-arm 18, an x-ray source 20, an x-ray receiving section 22, an optional calibration and tracking target 24 and optional radiation sensors 26. The optional calibration and tracking target 24 includes calibration markers 28 (see Figures 2a-2b), further discussed herein. A Carm controller 30 captures the x-ray images received at the receiving section 22 and stores the images for later use. The C-arm controller 30 may also control the rotation of the C-arm 18. For example, the C-arm 18 may move in the direction of arrow 32 or rotate about the long axis of the patient 14, allowing anterior or lateral views of the patient 14 to be imaged. Each of these movements involve rotation about a mechanical axis 34 of the C-arm 18. In this example, the long axis of the patient 14 is substantially in line with the mechanical axis 34 of the C-arm 18. This enables the C-arm 18 to be rotated relative to the patient 14, allowing images of the patient 14 to be taken from multiple directions or about multiple planes. An example of a fluoroscopic C-arm x-ray imaging device 16 is the "Series 9600 Mobile Digital Imaging System," from OEC Medical Systems, Inc., of Salt Lake City, Utah. Other exemplary fluoroscopes include biplane fluoroscopic systems, ceiling fluoroscopic systems, cath-lab fluoroscopic systems, fixed C-arm fluoroscopic systems, etc.

[0040] In operation, the imaging device 16 generates x-rays from the x-ray source 20 that propagate through the patient 14 and calibration and/or tracking target 24, into the x-ray receiving section 22. The receiving section 22 generates an image representing the intensities of the received x-rays. Typically, the receiving section 22 includes an image intensifier that first converts the x-rays to visible light and a charge coupled device (CCD) video camera that converts the visible light into digital images. Receiving section 22 may also be a digital de-

vice that converts x-rays directly to digital images, thus potentially avoiding distortion introduced by first converting to visible light. With this type of digital C-arm, which is generally a flat panel device, the calibration and/or tracking target 24 and the calibration process discussed below may be eliminated. Also, the calibration process may be eliminated for different types of medical procedures. Alternatively, the imaging device 16 may only take a single image with the calibration and tracking target 24 in place. Thereafter, the calibration and tracking target 24 may be removed from the line-of-sight of the imaging device 16.

[0041] Two dimensional fluoroscopic images taken by the imaging device 16 are captured and stored in the Carm controller 30. These images are forwarded from the C-arm controller 30 to a controller or work station 36 having the display 10 that may either include a single display 10 or a dual display 10 and a user interface 38. The work station 36 provides facilities for displaying on the display 10, saving, digitally manipulating, or printing a hard copy of the received images, as well as the five or six degree of freedom display. The user interface 38, which may be a keyboard, joy stick, mouse, touch pen, touch screen or other suitable device allows a physician or user to provide inputs to control the imaging device 16, via the C-arm controller 30, or adjust the display settings, such as safe zones of the display 10, further discussed herein. The work station 36 may also direct the C-arm controller 30 to adjust the rotational axis 34 of the C-arm 18 to obtain various two-dimensional images along different planes in order to generate representative two-dimensional and three-dimensional images. When the x-ray source 20 generates the x-rays that propagate to the x-ray receiving section 22, the radiation sensors 26 sense the presence of radiation, which is forwarded to the C-arm controller 30. to identify whether or not the imaging device 16 is actively imaging. This information is also transmitted to a coil array controller 48, further discussed herein. Alternatively, a person or physician may manually indicate when the imaging device 16 is actively imaging or this function can be built into the x-ray source 20, x-ray receiving section 22, or the control computer 30.

[0042] Fluoroscopic C-arm imaging devices 16 that do not include a digital receiving section 22 generally require the calibration and/or tracking target 24. This is because the raw images generated by the receiving section 22 tend to suffer from undesirable distortion caused by a number of factors, including inherent image distortion in the image intensifier and external electromagnetic fields. An empty undistorted or ideal image and an empty distorted image are shown in Figures 2a and 2b, respectively. The checkerboard shape, shown in Figure 2a, represents the ideal image 40 of the checkerboard arranged calibration markers 28. The image taken by the receiving section 22, however, can suffer from distortion, as illustrated by the distorted calibration marker image 42, shown in Figure 2b.

ed to three-dimensional volumetric images that can be

[0043] Intrinsic calibration, which is the process of correcting image distortion in a received image and establishing the projective transformation for that image, involves placing the calibration markers 28 in the path of the x-ray, where the calibration markers 28 are opaque or semi-opaque to the x-rays. The calibration markers 28 are rigidly arranged in pre-determined patterns in one or more planes in the path of the x-rays and are visible in the recorded images. Because the true relative position of the calibration markers 28 in the recorded images are known, the C-arm controller 30 or the work station or computer 36 is able to calculate an amount of distortion at each pixel in the image (where a pixel is a single point in the image). Accordingly, the computer or work station 36 can digitally compensate for the distortion in the image and generate a distortionfree or at least a distortion improved image 40 (see Figure 2a). A more detailed explanation of exemplary methods for performing intrinsic calibration are described in the references: B. Schuele, et al., "Correction of Image Intensifier Distortion for Three-Dimensional Reconstruction," presented at SPIE Medical Imaging, San Diego, California, 1995; G. Champleboux, et al., "Accurate Calibration of Cameras and Range Imaging Sensors: the NPBS Method," Proceedings of the IEEE International Conference on Robotics and Automation, Nice, France, May, 1992; and U.S. Patent No. 6,118,845, entitled "System And Methods For The Reduction And Elimination Of Image Artifacts In The Calibration Of X-Ray Imagers," issued September 12, 2000, the contents of which are each hereby incorporated by reference.

[0044] While the fluoroscopic C-arm imaging device 16 is shown in Figure 1, any other alternative imaging modality may also be used or an image-less based application may also be employed, as further discussed herein. For example, isocentric fluoroscopy, bi-plane fluoroscopy, ultrasound, computed tomography (CT), multi-slice computed tomography (MSCT), magnetic resonance imaging (MRI), high frequency ultrasound (HIFU), optical coherence tomography (OCT), intravascular ultrasound (IVUS), 2D, 3D or 4D ultrasound, -intraoperative CT, MRI, or O-arms having single or multi \*flat panels receivers that move about the ring to acquire 'fluoroscopic images may also be used to acquire preoperative or real-time images or image data of the patient 14. Image datasets from hybrid modalities, such as positron emission tomography (PET) combined with CT, or single photon emission computer tomography (SPECT) combined with CT, could also provide functional'image data superimposed onto anatomical data to be used to confidently reach target sights within the areas of interest. It should further be noted that the fluoroscopic C-arm imaging device 16, as shown in Figure 1, provides a virtual bi-plane image using a single-head C-arm fluoroscope 16 by simply rotating the C-arm 18 about at least two planes, which could be orthogonal planes to generate two-dimensional images that can be convert-

displayed on the six degree of freedom display 10. [0045] The navigation system 12 further includes an electromagnetic navigation or tracking system 44 that includes a transmitter coil array 46, the coil array controller 48, a navigation probe interface 50, an instrument 52 having an electromagnetic tracker and a dynamic reference frame 54. It should further be noted that the entire tracking system 44 or parts of the tracking system 44 may be incorporated into the imaging device 16, including the work station 36 and radiation sensors 26. Incorporating the tracking system 44 will provide an integrated imaging and tracking system. Any combination of these components may also be incorporated into the imaging system 16, which again can include a fluoroscopic C-arm imaging device or any other appropriate imaging device. Obviously, if an image-less procedure is performed, the navigation and tracking system 44 will

be a stand alone unit.

[0046] The transmitter coil array 46 is shown attached to the receiving section 22 of the C-arm 18. However, it should be noted that the transmitter coil array 46 may also be positioned at any other location as well, particularly if the imaging device 16 is not employed. For example, the transmitter coil array 46 may be positioned at the x-ray source 20, within the OR table 56 positioned below the patient 14, on siderails associated with the OR table 56, or positioned on the patient 14 in proximity to the region being navigated, such as by the patient's pelvic area. The transmitter coil array 46 includes a plurality of coils that are each operable to generate distinct electromagnetic fields into the navigation region of the patient 14, which is sometimes referred to as patient space. Representative electromagnetic systems are set forth in U.S. Patent No. 5,913,820, entitled "Position Location System," issued June 22, 1999 and U.S. Patent No. 5,592,939, entitled "Method and System for Navigating a Catheter Probe," issued January 14, 1997, each of which are hereby incorporated by reference.

[0047] The transmitter coil array 46 is controlled or driven by the coil array controller 48. The coil array controller 48 drives each coil in the transmitter coil array 46 in a time division multiplex or a frequency division multiplex manner. In this regard, each coil may be driven separately at a distinct time or all of the coils may be driven simultaneously with each being driven by a different frequency. Upon driving the coils in the transmitter coil array 46 with the coil array controller 48, electromagnetic fields are generated within the patient 14 in the area where the medical procedure is being performed, which is again sometimes referred to as patient space. The electromagnetic fields generated in the patient space induces currents in sensors 58 positioned in the instrument 52, further discussed herein. These induced signals from the instrument 52 are delivered to the navigation probe interface 50 and subsequently forwarded to the coil array controller 48. The navigation probe interface 50 provides all the necessary electrical isolation for the navigation system 12. The navigation probe interface 50 also includes amplifiers, filters and buffers required to directly interface with the sensors 58 in instrument 52. Alternatively, the instrument 52 may employ a wireless communications channel as opposed to being coupled directly to the navigation probe interface 50.

[0048] The instrument 52 is equipped with at least one, and may include multiple localization sensors 58. In this regard, the instrument 52 may include an orthogonal pair coil sensor 58 or a tri-axial coil sensor 58 or multiple single coil sensors 58 positioned about the instrument 52. Here again, the instrument 52 may be any type of medical instrument or implant. For example, the instrument may be a catheter that can be used to deploy a medical lead, be used for tissue ablation, or be used to deliver a pharmaceutical agent. The instrument 52 may also be an orthopedic instrument, used for an orthopedic procedure, such as reamers, impactors, cutting blocks, saw blades, drills, etc. The instrument 52 may also be any type of neurovascular instrument, cardiovascular instrument, soft tissue instrument, etc. Finally, the instrument 52 may be an implant that is tracked, as well as any other type of device positioned and located within the patient 14. These implants can include orthopedic implants, neurovascular implants, cardiovascular implants, soft tissue implants, or any other devices that are implanted into the patient 14. Particularly, implants that are formed from multiple components where the location and orientation of each component is dependent upon the location and orientation of the other component, such that each of these components can be tracked or navigated by the navigation and tracking system 44 to be displayed on the six degree of freedom display 10.

[0049] In an alternate embodiment, the electromagnetic sources or generators may be located within the instrument 52 and one or more receiver coils may be provided externally to the patient 14 forming a receiver coil array similar to the transmitter coil array 46. In this regard, the sensor coils 58 would generate electromagnetic fields, which would be received by the receiving coils in the receiving coil array similar to the transmitter coil array 46. Other types of localization or tracking may also be used with other types of navigation systems, which may include an emitter, which emits energy, such as light, sound, or electromagnetic radiation, and a receiver that detects the energy at a position away from the emitter. This change in energy, from the emitter to the receiver, is used to determine the location of the receiver relative to the emitter. These types of localization systems include conductive, active optical, passive optical, ultrasound, sonic, electromagnetic, etc. An additional representative alternative localization and tracking system is set forth in U.S. Patent No. 5,983,126, entitled "Catheter Location System and Method," issued November 9, 1999, which is hereby incorporated by reference. Alternatively, the localization system may be a hybrid system that includes components from various systems.

[0050] The dynamic reference frame 54 of the electromagnetic tracking system 44 is also coupled to the navigation probe interface 50 to forward the information to the coil array controller 48. The dynamic reference frame 54 is a small magnetic field detector or any other type of detector/transmitter that is designed to be fixed to the patient 14 adjacent to the region being navigated so that any movement of the patient 14 is detected as relative motion between the transmitter coil array 46 and the dynamic reference frame 54. This relative motion is forwarded to the coil array controller 48, which updates registration correlation and maintains accurate navigation, further discussed herein. The dynamic reference. frame 54 can be configured as a pair of orthogonally oriented coils, each having the same center or may be configured in any other non-coaxial coil configuration. The dynamic reference frame 54 may be affixed externally to the patient 14, adjacent to the region of navigation, such as the patient's spinal region, as shown in Figure 1 or on any other region of the patient. The dynamic reference frame 54 can be affixed to the patient's skin, by way of a stick-on adhesive patch. The dynamic reference frame 54 may also be removably attachable to fiducial markers 60 also positioned on the patient's body and further discussed herein.

[0051] Alternatively, the dynamic reference frame 54 may be internally attached, for example, to the spine or vertebrae of the patient using bone screws that are attached directly to the bone. This provides increased accuracy since this will track any motion of the bone. Moreover, multiple dynamic reference frames 54 may also be employed to track the position of two bones relative to a joint. For example, one dynamic reference frame 54 may be attached to a first vertebra. while a second dynamic reference frame 54 may be attached to a second vertebra. In this way, motion of the spine or vertebrae may be detected by the dual dynamic reference frames 54. An exemplary dynamic reference frame 54 and fiducial marker 60, is set forth in U.S. Patent No. 6,381,485, entitled "Registration of Human Anatomy Integrated for Electromagnetic Localization," issued April 30, 2002, which is hereby incorporated by reference.

[0052] Briefly, the navigation system 12 operates as follows. The navigation system 12 creates a translation map between all points in the radiological image generated from the imaging device 16 and the corresponding points in the patient's anatomy in patient space. After this map is established, whenever a tracked instrument 52 is used, the work station 36 in combination with the coil array controller 48 and the C-arm controller 30 uses the translation map to identify the corresponding point on the pre-acquired image, which is displayed on display 10. This identification is known as navigation or localization. An icon representing the localized point or instrument is shown on the display 10, along with five or six degrees of freedom indicia.

[0053] To enable navigation, the navigation system 12 will detect both the position of the patient's anatomy 14 and the position of the surgical instrument 52. Knowing the location of these two items allows the navigation system 12 to compute and display the position of the instrument 52 in relation to the patient 14. The tracking system 44 is employed to track the instrument 52 and the anatomy simultaneously. While the display 10 is configured to show the instrument with six degree of freedom accuracy.

[0054] The tracking system 44 essentially works by positioning the transmitter coil array 46 adjacent to the patient space to generate a low-energy magnetic field generally referred to as a navigation field. Because every point in the navigation field or patient space is associated with a unique field strength, the electromagnetic tracking system 44 can determine the position of the instrument 52 by measuring the field strength at the sensor 58-location. The dynamic reference frame 54 is fixed to the patient 14 to identify the location of the patient 14 in the navigation field. The electromagnetic tracking system 44 continuously recomputes the relative position of the dynamic reference frame 54 and the instrument 52 during localization and relates this spatial information to patient registration data to enable image guidance of the instrument 52 within the patient 14.

[0055] Patient registration is the process of determining how to correlate the position of the instrument 52 on the patient 14 to the position on the diagnostic, pre-acquired, or real-time images. To register the patient 14, the physician or user will select and store particular points from the pre-acquired images and then touch the corresponding points on the patient's anatomy with a pointer probe 62. The navigation system 12 analyzes the relationship between the two sets of points that are selected and computes a match, which correlates every point in the image data with its corresponding point on the patient's anatomy or the patient space. The points that are selected to perform registration are the fiducial arrays or landmarks 60. Again, the landmarks or fiducial points 60 are identifiable on the images and identifiable and accessible on the patient 14. The landmarks 60 can be artificial landmarks 60 that are positioned on the patient 14 or anatomical landmarks 60 that can be easily identified in the image data. The system 12 may also perform 2D to 3D registration by utilizing the acquired 2D images to register 3D volume images by use of contour algorithms, point algorithms, normalized mutual information, pattern intensity, or density comparison algorithms, as is known in the art.

[0056] In order to maintain registration accuracy, the navigation system 12 continuously tracks the position of the patient 14 during registration and navigation. This is necessary because the patient 14, dynamic reference frame 54, and transmitter coil array 46 may all move during the procedure, even when this movement is not desired. Therefore, if the navigation system 12 did not track the position of the patient 14 or area of the anato-

my, any patient movement after image acquisition would result in inaccurate navigation within that image. The dynamic reference frame 54 allows the electromagnetic tracking device 44 to register and track the anatomy. Because the dynamic reference frame 54 is rigidly fixed to the patient 14, any movement of the anatomy or the transmitter coil array 46 is detected as the relative motion between the transmitter coil array 46 and the dynamic reference frame 54. This relative motion is communicated to the coil array controller 48, via the navigation probe interface 50, which updates the registration correlation to thereby maintain accurate navigation.

[0057] It should also be understood that localization and registration data may be specific to multiple targets. For example, should a spinal procedure be conducted, each vertebra may be independently tracked and the corresponding image registered to each vertebra. In other words, each vertebra would have its own translation map between all points in the radiological image and the corresponding points in the patient's anatomy in patient space in order to provide a coordinate system for each vertebra being tracked. The tracking system 44 would track any motion in each vertebra by use of a tracking sensor 58 associated with each vertebra. In this way, dual displays 10 may be utilized, further discussed herein, where each display tracks a corresponding vertebra using its corresponding translation map and a surgical implant or instrument 52 may be registered to each vertebra and displayed on the display 10 further assisting an alignment of an implant relative to two articulating or movable bones. Moreover, each separate display in the dual display 10 may superimpose the other vertebra so that it is positioned adjacent to the tracked vertebra thereby adding a further level of information on the six degree of freedom display 10.

[0058] As an alternative to using the imaging system 16, in combination with the navigation and tracking system 44, the five or six degree of freedom alignment display 10 can be used in an imageless manner without the imaging system 16. In this regard, the navigation and tracking system 44 may only be employed and the probe 62 may be used to contact or engage various landmarks on the patient. These landmarks can be bony landmarks on the patient, such that upon contacting a number of landmarks for each bone, the workstation 36 can generate a three-dimensional model of the bones. This model is generated based upon the contacts and/or use of atlas maps. The workstation 36 may also generate a center axis of rotation for the joint or planes, based upon the probe contacts. Alternatively, the tracking sensor 58 may be placed on the patient's anatomy and the anatomy moved and correspondingly tracked by the tracking system 44. For example, placing a tracking sensor 58 on the femur and fixing the pelvis in place of a patient and rotating the leg while it is tracked with the tracking system 44 enables the work station 36 to generate a center of axis of the hip joint by use of kinematics and motion analysis algorithms, as is known in the art. If the

pelvis is not fixed, another tracking sensor 58 may be placed on the pelvis to identify the center of axis of the hip joint. If a tracking sensor 58 is placed on the femurand a tracking sensor 58 is placed on the tibia, upon moving this portion of the anatomy, a center of axis of the knee joint may be identified. Likewise, by placing a separate tracking sensor 58 on two adjacent vertebra and articulating the spine, the center of axis of the spinal region can also be identified. In this way, a target and/or model based on the center of the particular joint may be designated and identified on the six degree of freedom display 10. Movement of the instrument or implant 52 may then be tracked in relation to this target and/or model to properly align the instrument or implant 52 relative to the target and/or model.

[0059] Turning to Figures 3a and 3b, the method of employing the six degree of freedom display 10 is described in further detail. The method 64 begins by determining whether an image based medical procedure will be employed or an image-less medical procedure will be employed. If the image based procedure is being employed, the method proceeds along the first branch. In this regard, when an image based procedure will be utilized, the method begins at block 66 identifying the image tracking procedure. From block 66, the method proceeds to block 68 where images are generated by the imaging system 16. This imaging is performed at the area of interest of the patient 14 by any type of imaging device as previously discussed. Once images have been generated at block 68, the method proceeds to block 70 where calibration and registration is performed. In block 70, calibration of the imaging device 16 takes place using the calibration targets 28. Additionally, registration of the pre-acquired images from block 68 are registered to the patient space of the medical procedure utilizing the fiducial markers 60 and probe 62 as previously discussed. This registration registers the current patient space with the pre-acquired image, so that the instrument 52 or other devices may be tracked during the medical procedure and accurately superimposed over the pre-acquired images generated from the imaging device 16.

[0060] If an image-less medical procedure is selected, the method begins at block 72 identifying that an image-less based medical procedure will be performed. This method proceeds to either block 74 identifying a first way to generate image-less models or block 76 identifying a second way to generate image-less models. At block 74, the probe 62 is used to contact the body at various anatomical landmarks in the area of interest, such as a bone. For example, by touching the probe 62 to the pelvis, knee, and ankle, articulation planes can be defined using known algorithms and the center of each joint may also be defined. An example of this type of modeling is set forth in U.S. Patent No. 5,682,886, which is hereby incorporated by reference. Alternatively, multiple anatomical landmarks can be contacted with the probe 62 to generate a 3-D model with the more

points contacted, the more accurate the model depicted. [0061] Secondly, to generate a model at block 76, a tracking device is placed on the body and the body rotated about the joint. When this is done, the plane of rotation and joint center can be identified using known kinematic and/or motion analysis algorithms or using atlas maps or tables, as is known in the art. Once the area of interest has been probed, via block 74 or block 76, a model is generated at block 78. This model can be a 3-D surface rendered model, a 2-D model identifying articulating planes or a 3-D model identifying articulating planes and rotation, as well as the center of the joints. This enables the display 10 to use the joint centers or articulating planes as the target or trajectory, further discussed herein. 15

[0062] With each of the procedures 74 or 76, the procedure may be initially based on the use of atlas information or a 3-D model that is morphed, to be a patient specific model. In this regard, should the femur be the area of interest, an accurate representation of an ordinary femur may be selected from an atlas map, thereby providing an initial 2-D or 3-D model representing a typical anatomical femur. As with block 74, upon contacting numerous areas on the actual femur with the probe 62, the atlas model may be morphed into a patient specific 3-D model, with the more points contacted, the more accurate the morphed model. Patient specific information may also be acquired using an ultrasound probe to again identify the shape of the patient's natural femur in order to morph the atlas model. A fluoroscopic image of the region may also be used to morph the patient's femur with the atlas model to provide a patient specific morphed model. Proceeding under block 76 and assuming that the area of interest is the hip joint, an atlas model of the femur and pelvis may be the initial starting point. Upon rotating and moving the femur relative to the pelvis, a patient specific morphed model may be created to generate accurate joint centers and axes of motion again using known kinematics and/or motion analysis algorithms

[0063] Once the image data is calibrated and registered at block 70 or the model is generated at block 78, the method proceeds to block 80. At block 80, the specific type of coordinate system is selected, which will be displayed by indicia on the six degree of freedom display 10. The coordinate systems can be a Cartesian coordinate system, a spherical coordinate system, or a polar coordinate system. By way of example, the Cartesian coordinate system will be selected. The Cartesian coordinate system will include the X, Y, and Z axes, and X rotation, Y rotation, and Z rotation about its respective axes.

[0064] With reference to Figure 5, the six degree of freedom display 10 is shown in further detail employing the Cartesian coordinate system. In this regard, the X axis 82 and the Y axis 84 are shown positioned on the display 10. The Z axis 86 extends out from the display 10 and is shown in the upper left corner. Rotation about

the X axis is shown by bar graph 88 to the right of the display 10 and rotation about the Y axis is shown by bar graph 90 positioned at the bottom of the display 10. Rotation about the Z axis is shown with the arcuate bar graph 92 oriented about the X and Y axes 82 and 84. Each axis, as well as the rotation axes identified by the bar graphs may be color coded to identify safe zones or regions for the item being tracked or navigated. In this regard, the safe zones can be defined as ranges around the planned trajectory path or target where the safe zones are determined by manufactured determined parameters, user determined parameters or patient specific parameter, further discussed herein.

[0065] Arrow indicator 94 identifies the degree of ro-

tation about the X axis 82. Arrow indicator 96 shows the amount of rotation about the Y axis 84. Arrow 98 identifies the rotation about the Z axis, while arrow 100 identifies the depth being tracked along the Z axis 86. The \* origin 102 may be set to be the desired target position or trajectory path. The crosshairs 104 represents the tip is following the instrument 52 being tracked, while the circle 106 represents the hind area of the instrument 52 being tracked. With the understanding that the instrument 52 can be any type of medical device or implant. Also, if five degree of freedom information is provided, one of the indicia 82, 84, 86, 88, 90, and 92 will be removed. [0066] Once the coordinate system is selected at block 80, the method proceeds to block 108 where the target or trajectory is selected. The target or trajectory selected at block 108 is typically positioned at the origin 102 on the display 10. In this way, the object being tracked or aligned may be tracked and aligned about the origin 102. Alternatively, the target may be identified at any coordinate within the display 10 or multiple targets may also be identified within the single display 10. An indicia of the target may also be positioned on the display 10. The target is selected based upon the desired area to position the instrument 52 and can be selected from the pre-acquired images or from the 3-D model. Once selected, this target is correlated to the display 10 and generally positioned at the origin 102.

\* [0067] Once the target/trajectory is selected at block 108, such as the origin 102, the method proceeds to degree of freedom. Referring again to Figure 5, the safe zones 112 are identified for each degree of freedom by color coding. For example, the safe zone 112 for the X axis is between -2.5 and +2.5. The safe zone 112 for rotation about the X axis is between -5° and +5° of ro-\*\*\* tation about the X axis. The user can simply guide the instrument 52 using the cross hairs 104 and circle 106 to align the instrument 52 within these designated safe \*\*\*\*Zones 112. Again, these safe zones 112 may be determined by manufacture specifications, such as tolerance range of the instruments or positions for implants. The safe zones 112 may also be determined based on the patient, the surgeon conducting the procedure, or any other factors to assist a surgeon in navigating the instrument 52 through the patient 14. These safe zones 112 may also be identified via an audible signal or tone or a varying tone. The safe zones 112 may also be identified by any other convenient manner to be set out on the display 10.

[0068] Once the safe zones 112 are identified for each degree of freedom in block 110, the method proceeds to block 114 where the target trajectory in the selected coordinate system is displayed with the safe zones 112, as shown in Figure 5. At block 116, if an image based method is being conducted, a decision whether to superimpose the image over the target/trajectory is made. Alternatively, the image may be placed adjacent to the target trajectory display, as is shown in Figure 6, and further discussed herein. Should the image-less based. medical procedure be conducted, at block 118, a determination is made whether to superimpose the model that was generated at block 78. Here again, this model may be superimposed over the target/trajectory display on display 10 or may be positioned adjacent to the target/trajectory display in a split screen or on a separate display.

[0069] Once the target/trajectory 102 is displayed along with the safe zones 112 in the proper coordinate system, as shown in Figure 5, the method proceeds to block 120 where the first implant/instrument 52 is tracked with the navigation system 44. With the implant/ instrument 52 being tracked at block 120, the method proceeds to block 122 wherein indicia representing the implant/instrument 52 is displayed on the display 10, with either five or six degrees of freedom information. Here again, referring to Figure 5, the indicia representing the implant/instrument 52 is the crosshairs 104 and the circle 106 designating the tip and hind, respectively. The tip 104 and hind 106 is represented in relation to the target/trajectory path 102 in six degrees of freedom. This six degrees of freedom include the X and Y locations, as well as the depth Z of the implant/instrument 52 displayed. In addition, the rotation about each of the axes is also provided. This rotation can be helpful in many applications, including orthopedic, where rotation specific components need to be positioned relative to one another. For example, in a spinal application, alignment of a pedicle screw in relation to a spinal rod, would require information regarding the rotation of the screw relative to the rod. In cardiac procedures, this may be useful where ablation is necessary on a certain side of an artery and the ablation electrode is only positioned on one side of the catheter. In this situation, rotation of the catheter relative to the target in the artery is critical. In a neuro procedure, a biopsy needle may only have a biopsy port positioned upon one portion of the circumference of the needle, thereby requiring the rotation of the biopsy needle to be known in order to provide the proper capture of the relevant biopsy sample. Without this display, this information would not be available. [0070] With the indicia of the implant/instrument 52

being displayed, the implant/instrument 52 is aligned or

9

fixed with the target/trajectory 102 at block 124. In this regard, the tip 104 and the hind 106 are aligned and fixed relative to the target/trajectory 102 at the origin and the rotational orientation is also aligned to the desired position. Again, the target/trajectory 102 may not be positioned at the origin and can be positioned anywhere within the coordinate system if desired. As shown in Figure 5, the tip 104 of the implant/instrument 52 is shown aligned with the target 102, while the hind 106 is slightly offset from the target/trajectory 102. Once the implant/instrument 52 is aligned and fixed relative to the target/trajectory 102, the method proceeds to block 126.

[0071] At block 126, a determination is made as to whether there is a second implant/instrument 52 to be tracked. If there is not a second implant/instrument 52 to be tracked, the method ends at block 128. Should there be a second implant/instrument 52 to track, such as a corresponding implant component that articulates with the first implant, the method proceeds to block 130. At block 130, a second target/trajectory 102 is selected, which is based upon the alignment or fixation of the first implant/instrument 52 relative to the first target/trajectory 102. In this regard, if the surgeon is not able to position the first implant/instrument 52 at the desired target/trajectory 102, this offset from the target/trajectory 102 may affect the second implant, which possibly articulates or mates with the first implant. If this is the case, the second target/trajectory 102 will need to take into consideration this offset in order to provide proper articulation and alignment of the first implant component with the second implant component.

[0072] With minimally invasive types of procedures, the implant may also have numerous components with each component articulating or mating with another component, thereby requiring tracking of each component as it is implanted during the minimally invasive procedure. This second target/trajectory 102 may be displayed on a separate display 10 (see Figure 1), positioned via a split screen of a single display 10 or may be superimposed upon the existing display that displays the first target 102 and implant position. In this way, orientation and placement of both the first and second implants, which are dependent upon one another can be shown in the display 10 providing the surgeon the opportunity to adjust either position of either implant intraoperatively before the implants are permanently affixed to the patient 14. These types of implants include knee implants, hip implants, shoulder implants, spinal implants, or any other type of implant, which has a bearing surface and an articulating surface or any type of implant having multiple mating and connecting components.

[0073] Once the second target/trajectory 102 has been selected at block 130, the method proceeds to block 132. At block 132, the safe zones 112 for each degree of freedom is selected for the second implant/instrument 52 similar to the way the first set of safe zones 112 were selected for the first implant/instrument 52. Once the second safe zones 112 are selected, the

method proceeds to block 134. At block 134, the display 10 displays the second target/trajectory 102 in the same coordinate system with the second safe zones 112. Here again, at block 136, if it is an image based medical procedure; the pre-acquired image may be superimposed on to the target/trajectory 102. Alternatively, this image can be positioned adjacent the target screen in a split screen configuration (see Figure 6 and 7). If the method is proceeding as an image-less type medical procedure, at block 138, decision is made whether to superimpose the generated model from block 78. Once the target/trajectory 102 is in the proper coordinate system with the safe zone 112 are displayed at display 10, the surgical implant/instrument 52 is tracked at block 140. Here again, the second implant/instrument 52 can be tracked on a separate display 10 or be tracked on the same display as the first implant/instrument 52.

[0074] Alternatively, separate displays 10 may be used where information is linked between the displays showing the second implant/instrument 52 in relation to the first implant/instrument 52. With the second implant/ instrument 52 being tracked at block 140, the second implant/instrument 52 is displayed in relation to the second target/trajectory 102 in five or six degrees of freedom at block 142. Again, this may be a separate display 10, a split screen display 10 with both the first target/ trajectory 102 and the second target/trajectory 102 or the same display 10 displaying both targets/trajectories 102. While the second implant/instrument 52 is being displayed, the second implant/instrument 52 is aligned and fixed at the second target/trajectory 102 at block 144. Once the second implant/instrument 52 is fixed at block 144, the method proceeds to block 146.

[0075] At block 146, a determination is made whether the alignment or fixation of the first and second implants/ instruments 52 are correct. In this regard, with two separate displays 10 linked or with a single display 10, showing both targets/trajectories 102, a surgeon can determine whether each implant/instrument 52 is within its desired safe zones 112 and, therefore, optimally positioned for proper articulation. Here again, these safe zones 112 may be color coded for the different safe zones provided. If both implants are positioned and fixed at the proper targets, the method ends at block 148. If one or both of the implants are not properly positioned, adjustment of the first or second target/trajectory 102 is performed at block 150. Once either or both targets are adjusted, realignment of the first and/or second implants/instruments 52 are performed at block 152. Here again, since multiple component implants are dependent upon one another with respect to their position and orientation, alignment and adjustments of the targets/ trajectories 102 may be performed several times until the optimum placement for each is performed at repeat block 154. Thereafter, the method terminates at end block 156.

[0076] While the above-identified procedure is discussed in relation to an orthopedic medical procedure

in which an implant having multiple implant components is implanted within a patient using the six degree of freedom display 10, it should be noted that the six degree of freedom display 10 may be used to track other medical devices as well. For example, as was briefly discussed, an ablation catheter generally has an electrode positioned only on one angular portion of its circumference. Likewise, the wall of an artery typically has a larger plaque build-up on one side. Therefore, it is desirable to align that ablation electrode with the proper side of the artery wall during the procedure. With the six degree of freedom display 10, the surgeon can easily identify the location, depth and angular rotation of the catheter relative to the artery wall. Other types of procedures may require the medical instrument or probe to be properly oriented and located within the patient, such as identifying and tracking tumors, soft tissue, etc. By knowing and displaying the six degree of freedom movement of the medical device on the display 10, the medical procedure is optimized.

[0077] It should also be pointed out that the method discussed above requires that the implant/instrument 52 have a tracking sensor associated therewith in order to identify the location of the tracked device in six degrees of freedom and display it on the display 10. The tracking sensors may be attached directly to implants, attached to the instruments that engage the implants or attach to members extending out from the implants. These tracking sensors again may be electromagnetic tracking sensors, optical tracking sensors, acoustic tracking sensors, etc. Examples of various targets, which may or may not be superimposed on the display again include orthopedic targets, spinal targets, cardiovascular targets, neurovascular targets, soft tissue targets, etc. Specific examples include again the location of the plaque on a wall of an artery, the center of an articulating joint being replaced, the center of the implant placement, etc. By displaying two targets, either on separate displays or on the same display, the surgeon can dynamically plan and trial implant placements by moving one component of the implant to see where the other articulating component of the implant should be positioned. In this way, the surgeon can trial the implant confirming its placement and orientation, via the display 10 before the implant is permanently affixed to the patient 14.

[0078] In a spinal procedure, two adjacent vertebra bodies can be tracked and displayed on two separate displays. In this way, if a single jig, such as a cutting jig is used to cut both the surface of the first vertebra and the surface of the second vertebra, orientation of the jig may be displayed on each separate display in relation to the corresponding vertebra being acted upon, thereby enabling simultaneous tracking of the two planes being resected for each separate vertebra on a dual display system. Additionally, each vertebra may be displayed on each of the dual displays so that the vertebra being tracked is shown with the adjacent vertebra superim-

posed adjacent thereto. Once the vertebra bodies are prepared, the implant is typically placed between each vertebra on the prepared site. Other ways of preparing this site is by using drills, reamers, burrs, trephines or any other appropriate cutting or milling device.

[0079] Briefly, the method, as shown in Figures 3a and 3b, demonstrates that the display 10 illustrated both the position and orientation of an object with respect to a desired position and orientation with six degrees of freedom accuracy. The display 10 may be automatically updated in real-time using the navigation system 44 to report the orientation of the tracked device. The user may also adjust the display 10 in order to control a device's orientation. The display 10 again consists of three rotational indicators (RX, RY, RZ) and three translational indicators or indicia (TX, TY, TZ). Each indicator shows both visual and quantitative information about the orientation of the device. Each indicator also displays a predetermined safe zone 112 and application-specific label for each degree of freedom. As noted, it may also be relevant to overlay the display 10 over anatomical image data from the imaging device 16. When working with 3-D image data sets, the anatomy normal to the tip 104 of the positioned device can provide the user with additional positional information.

[0080] Tones, labels, colors, shading, overlaying with image data can all be modified and incorporated into the display 10. The current display 10 is also shown as a Cartesian coordinate based display, but again could be based on a polar based display or a spherical based display and a quick switch between both can be supplied or simultaneously displayed. The display can also be configured by the user to hide parameters, location, size, colors, labels, etc.

[0081] Some medical applications that may be commonly displayed and linked to the display 10 are: 1) reaming of an acetabular cup with major focus upon RY and RZ, 2) length of leg during hip and knee procedures focused upon TZ and RZ, 3) biopsies and ablations focused upon RX, RY, and RZ for direction of the therapy device, and 4) catheters with side ports for sensing information or delivery of devices, therapies, drugs, stem cells, etc. focused upon six degree of freedom information.

[0082] Referring now to Figures 4a - 4e, a medical procedure employing a six degree of freedom alignment display 10 is shown in further detail. In this example, an orthopedic medical procedure replacing the hip joint is illustrated. During this procedure, various instruments 52, as well as the implants 52 are tracked and aligned using the six degree of freedom display 10. Referring specifically to Figure 4a, a femur 160 having a femoral head 162 is illustrated, along with a pelvis 164 having an acetabulum 166. Assuming that the medical procedure being performed is an image based system, this area of interest will be imaged by the imaging device 16. Here again, the dynamic reference frame 54 may be attached to the femur 154 or the pelvis 164 or two dynamic

frames 54 may be attached, one to each bone to provide additional accuracy during the medical procedure. With the head 162 dislocated from the acetabulum 166, a center of articulation of the acetabulum 166 is identified as the target 168, shown in Figure 6.

[0083] In this regard, Figure 6 illustrates the display 10 configured as a split screen with the right identifying the six degree of freedom display and the left illustrating the pre-acquired image with the center of articulation 168 being the intersection of the X, Y, and Z axes. As illustrated in Figure 4a, a reamer 170 having a tracking sensor 58 is shown reaming the acetabulum 166. The tracking system 44 is able to accurately identify the navigation of the tip and hind of the reamer 170. As illustrated in Figure 6, in the right half of the split screen, one can observe that the tip represented by the crosshairs 172 is properly positioned along the X and Y coordinates and within the corresponding safe zones 112, however, the hind portion of the instrument 170, as identified by the circle 174, is angularly offset from the target 168 at the origin. The surgeon can then angularly adjust the hind portion 174 of the instrument 170 until the hind portion 174 is shown in the display 10 as positioned over the crosshairs 172, thereby assuring proper alignment of the reaming device 170 for subsequent proper placement of the acetabular cup implant. By tracking the reamer 170, the surgeon can be relatively confident that an acetabular cup implant will be properly positioned before the implant is even impacted into the acetabulum 166.

Turning to Figure 4b, an acetabular cup 178 is [0084] shown being impacted into the reamed acetabulum 166, via the tracked guide tool 180 with an impactor 182. The guide tool 180 has a distal end, which is nestingly received within the acetabular cup 178. Thus, by tracking the instrument 180, via tracking sensor 58, orientation of the acetabular cup 178 may be displayed on the display 10 in six degrees of freedom. In this way, before the acetabular cup 178 is impacted into the acetabulum 166, the surgeon can view on the display 10 whether the acetabular cup 178 is properly positioned at the proper angular orientation, as shown in Figure 7, the impactor 180 is shown superimposed over an image generated by the imaging device 16. In this way, the proper orientation, including abduction and anteversion is achieved before the acetabular cup 178 is permanently implanted.

[0085] Once the acetabular cup 178 has been impacted, the femoral head 162 is resected along a plane 184 by use of a cutting guide 186, having the tracking sensor 58 and a saw blade 188. By using the center of the femoral head 162 as the second target, the cutting plane 184 may be properly defined to provide proper articulation with the acetabular cup 178 before a hip stem is implanted in the femur 160. Here again, the second target is dependent upon the first target. Thus, if the acetabular cup 178 was implanted somewhat offset from its target. Hie second target may be properly com-

pensated to accommodate for this offset by use of the display 10. In this regard, a second display illustrating the target for the cutting plane 184 may be provided.

[0086] Once the femoral head 162 of the femur 160 has been resected, as shown in Figure 4d, a reamer 190 is employed to ream out the intramedullary canal 192 of the femur 160. In order to provide proper angular orientation of the reamer 190, as well as the depth, a subsequent target can be defined and identified on the display 10 and tracked by use of the tracking sensor 58. This target may be displayed separately or in combination with the previously acquired targets. By insuring the proper angle of the reamer 190 relative to the longitudinal axis of the femur 160 is tracked and displayed on display 10, the surgeon can be provided a higher level of confidence that the hip stem will be properly positioned within the intramedullary canal 192.

[0087] Once the intramedullary canal 192 has been reamed by the reamer 190, a hip stem 194 is impacted with an impactor 196 into the intramedullary canal 192. By targeting the acetabular cup location, along with the resection plane 184 and the reaming axis of the reamer 190, upon positioning the hip stem 194, within the femur 160, proper articulation and range of motion between the acetabular cup 178 and the hip stem 194 is achieved without time consuming trialing as is conducted in conventional orthopedic procedures. Thus, by providing the safe zones 112 in relation to the hip stem 194 size, proper articulation with the acetabular cup 178 is achieved. Here again, while an example of an orthopedic hip replacement is set out, the six degree of freedom display 10 may be utilized with any type of medical procedure requiring visualization of a medical device with six degree freedom information.

[0088] The six degree of freedom display 10 enables implants, devices and therapies that have a specific orientation relative to the patient anatomy 14 to be properly positioned by use of the display 10. As was noted, it is difficult to visualize the correct placement of devices that require five or six degree of freedom alignment. Also, the orientation of multiple-segment implants, devices, or therapies in five and six degrees of freedom so that they are placed or activated in the correct orientation to one another is achieved with the display 10. Since the location and orientation is dependent upon one another to be effective, by having the proper orientation, improved life of the implants, the proper degrees of motion, and patient outcome is enhanced. Also, the six degree of freedom display 10 may be used as a user input mechanism by way of keyboard 38 for controlling each degree of freedom of a surgical robotic device. In this regard, the user can input controls with the joystick, touch screen or keyboard 38 to control a robotic device. These devices also include drill guide holders, drill holders, mechanically adjusted or line devices, such as orthopedic cutting blocks, or can be used to control and drive the alignment of the imaging system 16, or any other type of imaging system.

Since multiple implants and therapies, or multi-[0089] segment/compartment implants require multiple alignments, the display 10 may include a stereo display or two displays 10. These displays may or may not be linked, depending on the certain procedure. The target point/location (translation and orientation of each implant component is dependent upon the other implant placement or location). Therefore, the adjustment or dynamic targeting of the dependent implant needs to be input to the dependent implant and visually displayed. Again, this can be done by two separate displays or by superimposing multiple targets on a single display. Many implants such as spinal disc implants, total knee and total hip replacements repair patient anatomy 14 by replacing the anatomy (bone, etc.) and restoring the patient 14 to the original biomechanics, size and kinematics. The benefit of the six degree of freedom alignment display 10 is that original patient data, such as the images can be entered, manually or collectively, via the imaging device 16 or image-less system used for placement of the implant. Again; manually, the user can enter data, overlay templates, or collect data, via the imaging system 16. An example, as discussed herein of an application is the alignment of a femoral neck of a hip implant in the previous patient alignment. The previous patient alignment can be acquired by landmarking the patient femoral head by using biomechanics to determine the center and alignment of the current line and angle of the femoral head. This information can be used as the target on the display 10 in order to properly align the implant replacing the femoral head.

[0090] The six degree of freedom display 10 also provides orientation guidance on a single display. Separate visual and quantitative read-outs for each degree of freedom is also displayed on the display 10. Visual representations or indicia of procedure-specific accepted values (i.e., a "safe zone 112") for each degree of freedom is also clearly displayed on the display 10. These safe zones 112 are displayed as specifics or ranges for the user to align or place within. The procedure specific accepted values for the safe zones 112 can be manufacture determined, user determined, patient specific (calculated) or determined from algorithms (finite ele-2 ment analysis, kinematics, etc. atlas or tables). It can also be fixed or configurable. Safe zones 112 can also be defined as ranges around a planned trajectory path or the specific trajectory path itself (range zero). The trajectory paths are input as selected points by the user or paths defined from the patient image data (segmented vascular structure, calculated centers of bone/joints, anatomical path calculated by known computed methods, etc.).

[0091] Turning now to Figures 8a-8g, another medical procedure that may employ the six degree of freedom alignment display 10 is shown in further detail, along with Figure 9 illustrating the use of the display 10 during this medical procedure. In this example, a spinal medical procedure that implants a cervical disc implant be-

tween two vertebrae is illustrated. During this procedure, various instruments 52, as well as the implant 52 are tracked and aligned using the six degree of freedom display 10. Also, the bony structures during the procedure are also tracked.

[0092] Referring specifically to Figure 8a, a first vertebra or vertebral body 200 is shown positioned adjacent to a second vertebra or vertebral body 202 in the cervical area of the spine. Assuming that the medical procedure is being performed in an image based system, this area of interest would be imaged by the imaging device 16. Again, a dynamic reference frame 54 may be attached to the first vertebra 200 and a second dynamic reference frame 54 may be attached to the second vertebra 202. These dynamic reference frames 54may also be combined with tracking sensors 58, which are shown attached to the vertebral bodies 200 and 202. A center of articulation of the vertebra 200 and a center of articulation of a vertebra 202 may be identified as the targets 168 on the dual display illustrated on Figure 9. In this way, by utilizing the center of articulation of each vertebral body with respect to each other as the targets 168, tracking of the instruments 52 used during the procedure, as well as the implant 52 with respect to these articulation centers may be achieved. This center of articulation or instantaneous center of rotation is identified as the "X" along axis Y. A plane or axis X is shown perpendicular to the longitudinal or spinal axis Y. This axis is where the implant, as well as milling should be performed or centered around.

[0093] Referring to Figure 8b, a cam distracter instrument 204 is shown distracting the vertebra 200 relative to the vertebra 202. The cam distracter 204 may be tracked, via another tracking sensor 58 affixed to the cam distracter 204. In this way, the six degree of freedom display 10 illustrated in Figure 9 can illustrate a location of the cam distracter 204 relative to the center of each vertebra 200 and 202 independently on the display. Since the instrument 204 is rigid, by locating the tracking sensor 58 on the instrument 204, the distal end of the instrument 204 is known and may be illustrated on the display 10 using crosshairs 104 and circle 106 to represent the tip and hind, respectively.

[0094] Once each vertebrae 200 and 202 have been distracted by the cam distracter 204, a sagittal wedge 206 also having a tracking sensor 58 is utilized and shown in Figure 8c. The sagittal wedge 206 is used to center each vertebrae 200 and 202, along the sagittal plane and again may be tracked and displayed with six degree of freedom on the display 10, as illustrated in Figure 9. In this regard, the surgeon can confirm both visually and via the display 10 that the sagittal wedge 206 is centered on the sagittal plane between the vertebrae 200 and 202, as well as obtain the proper depth, via the Z axis display 86 on the display 10, illustrated in Figure 9.

[0095] Once the sagittal centering has been achieved with the sagittal wedge 206, the medical procedure pro-

ceeds to burring as shown in Figure 8d. In this regard, a burr 208 attached to a burring hand piece 210, also having a tracking sensor 58, is used to burr an area between the first vertebra 200 and the second vertebra 202. Here again, the orientation of the burr 208 relative to each vertebrae 200 and 202 may be displayed on the display 10 with six degree of freedom information. Therefore, burring along the proper X and Y plane, as well as the proper depth may be visually displayed with the appropriate indicia, as illustrated in Figure 9. Rotational information about the corresponding X, Y and Z axes is also displayed. By burring within the safe zones 112 using the information regarding the surgical implant 52 as the safe zones 112, the surgeon can be assured to perform the proper burring between the vertebrae 200 and 202 to insure a proper oriented fit for the surgical implant 52. By tracking the burr 208 with six degrees of freedom information, the mounting anchors 212 for the hand piece 210 are optional and may not be required. Additionally, each single display in the dual display 10, as shown in Figure 9, may also superimpose an image of each vertebrae 200 and 202 relative to one another on the display with each display having its coordinate system referenced to one of the vertebrae. The resulting milled vertebrae 200 and 202 are shown in Figure 8e with a ring portion 214 milled to receive the spinal implant 52.

[0096] Referring to Figures 8f and 8g, the spinal implant 52 is shown being implanted between the vertebrae 200 and 202 using an implant inserter 216 that is also tracked by tracking sensor 58. The spinal implant 52 may be any type of cervical or other spinal disc implant for any other area of the spine. For example, the spinal implant may be the spinal implant disclosed in U. S. Patent No. 5,674,296, entitled "Human Spinal Disc Prosthesis," issued October 7, 1997, U.S. Patent No. 5,865,846, entitled "Human Spinal Disc Prosthesis," issued February 2, 1999, also known as the Bryan Cervical Disc System, offered by Medtronic Sofamor Danek of Minneapolis, Minnesota or the Prestige Cervical Disc System, also offered by Medtronic Sofamor Danek, or any other spinal disc implant, all of which are hereby incorporated by reference. By tracking the implant inserter 216 relative to the vertebrae 200 and 202, proper orientation of the spinal implant 52, as well as rotational orientation about the Z axis can be clearly displayed on the six degree of freedom display 10, as shown in Figure 9. Rotation about the Z axis is used to make sure that the flanges 218 of the implant 52 are properly oriented and centered along the sagittal plane, as shown in Figure 8g. Again, by using the display 10, as illustrated in Figure 9, the anchors 220 are optional since orientation... of the implant 52 can be tracked continuously as it is inserted between the vertebrae 200 and 202. Here again, this eliminates the need for forming holes in the vertebrae 200 and 202. It should further be noted that the implant 52 illustrated in these figures is merely an exemplary type of spinal implant and any known spinal

implants may also be similarly tracked. For example, another common type of spinal implant is formed from a two-piece unit that includes a ball and cup articulating structure that may likewise be independently tracked to assure proper fit and placement.

[0097] Here again, the six degree of freedom display 10, which is illustrated as a split or dual display 10 in Figure 9 assists a surgeon in implanting a spinal implant 52 in order to achieve proper fixation and orientation of the implant 52, relative to two movable vertebrae 200 and 202. By tracking each vertebra 200 and 202 independently, and tracking its resection, should one vertebra be resected off-plane due to anatomical anomalies, adjustment of the plane at the adjacent vertebra may be achieved in order to still provide a proper fit for the spinal implant 52. In this way, each vertebrae 200 and 202 can be independently monitored, so that if one is off axis, the other can be manipulated accordingly to account for this adjustment. Additionally, by monitoring the entire process having six degree of freedom information, via display 10, further accuracy was achieved, thereby providing increased range of motion for the patient after implantation of the implant 52.

[0098] By use of the six degree of freedom display, for the various types of medical procedures, improved results can be achieved by providing the surgeon with the necessary information required. In regard to surgical implants, the range of motion may be increased while reducing impingement of two-part articulating or fixed implants. This also enables maximum force transfer between the implant and the body. With therapy delivery procedures, by knowing the location of the catheter delivery tube and the specific port orientation, accurately aiming at the site is enabled to provide maximum delivery of the therapy at the correct site. This procedure also enhances and enables better results when using an ablation catheter by again knowing the rotational orientation of the ablation catheter and the ablation electrode relative to the area in the wall of the artery that requires ablation. Finally, by knowing the rotational orientation of a ablation or biopsy catheter, this type of catheter may be easily directed and aligned to tumors, stem cells, or other desired sites in an easy and efficient manner.

[0099] Turning to Figure 10, a method 230 for post-operative adjustment or tuning of implants, such as a spinal implant, according to the teachings of the present invention is disclosed. The method 230 also includes pre-operative planning, implanting, as well as the post-operative exam procedure. In this regard, the method 230 begins at block 232 where pre-operative planning of the medical procedure begins. The pre-operative planning proceeds from block 232 to either block 234 if an image based pre-operative plan is conducted or block 236 if both an image and sensing pre-operative plan is being conducted, the method proceeds to block 238 where pre-operative image data is acquired. The pre-operative images may be captured from a four-dimen-

sional CT scan, which provides for capturing images over a specific time frame. In this regard, if the pre-operative planning is for implantation of a cervical disc, the patient may be asked to move his or her neck in different manners to capture the image data over time. Alternatively, any other type of imaging device 16 may be employed to either simply gather static image data or image data over time. The captured image data may also be used in conjunction with the electromagnetic tracking system 44, as discussed herein. Another example of pre-operative planning using a tracking system is disclosed in U.S. Patent No. 6,470,207, entitled "Navigation Guidance Via Computer-Assisted Fluoroscopic Imaging," issued October 22, 2002, which is hereby incorporated by reference. Other types of pre-operative planning using a tracking system may also be employed. This image data is then analyzed at the analysis data block 240, further discussed herein.

[0100] Should the pre-operative planning proceed to block 236, which employs the image and sense-based pre-operative planning, this procedure will capture image data and sense parameters at block 242. In this regard, the captured image data may be the same image data that is captured at block 238. In addition to the captured image data, various parameters in the area of interest may also be sensed during the pre-operative planning state. In this regard, probes or sensors maybe placed in the area of interest to sense the parameters, such as temperature, pressure, strain, and force motions. For example, in a cervical disc implant, sensors may be positioned between adjacent vertebrae of interest to measure temperature in certain areas, which may indicate friction or impingement. Likewise, strain gauges may be positioned to measure forces to identify areas having unacceptably high forces between the vertebrae. Again, this data is then analyzed at block 240.

[0101] At block 240, the data from either the image based or the image sense based pre-operative planning is analyzed. Should the data only include image data from block 238, this image data may be used to identify areas of interest, the patient size, and be used to assist in preparing the surgical plan. By viewing this data, such as 4D data, which is essentially 3D data over time or static image data, certain abnormal or irregular movements in the area of interest may be identified. These areas may be identified by visual examination, by performing finite element analysis or other known motion analysis to create a 3D model of the captured image. The finite element analysis may include calculating the instantaneous center of rotation "x" or make this determination from the image data itself. The overall shape of the spine may also be analyzed via the image data to identify and determine various force vectors on the discs of interest by analyzing the entire spine, the curvature of the spine and the articulation area of the angle of the spine relative to the ground. This information may be used to find force vectors and loading on the various regions of the vertebrae of interest. Should the sensed parameters also be used, or alternatively only be used, these sensor readings, which can be measured statically or actively while the patient is moving are utilized to again identify points of interest or potential abnormal activities by sensing parameters, such as temperature, pressure, stress, and strain in the area of interest.

[0102] Once the data has been analyzed at block 240, the procedure proceeds to block 243, where the implant and the type of procedure is selected. The implant is selected, based on the various abnormalities identified in order to enable the surgeon to resolve the noted abnormalities. The implant is selected based on various parameters, such as material selection, performance characteristics, stiffness, style or implant type and sizing. Once the type of implant has been selected, sizing of the implant may also be pre-operatively performed, based on the data captured and analyzed at block 240. Sizing may be performed using known sizing templates, which provides the surgeon with a visual means of correlating the size of the implant to the area of interest. Alternatively, various sized templates automated in software may also be included and stored within the work station 36 and superimposed in the area of interest to provide a visual indication of the sized implant to select. In addition to selecting the type and size of the implant, the type of procedure to position the implant may be determined pre-operatively.

[0103] Once the size and type of implant is selected, as well as the type of procedure, the procedure proceeds to block 244. At block 244, the selected implant is implanted generally under surgical guidance in the area of interest. For example, a cervical disc implant may be implanted, as illustrated in Figures 8a-8g. However, any other type of implant procedure may also be performed to position the selected implant, which may include a non-surgically guided procedure. Other exemplary types of surgically guided procedures are set out in U.S. Patent No. 6,470,207, entitled "Navigation Guidance Via Computer-Assisted Fluoroscopic Imaging," issued October 22, 2002; U.S. Patent No. 6,434,415, entitled "System For Use in Displaying Images Of A Body Part," issued August 13, 2002; and U.S. Patent No. 5,592,939, entitled "Method and System for Navigating a Catheter Probe," issued January 14, 1997, all of which are hereby incorporated by reference.

[0104] After implantation, there is a recovery period, exemplified by block 246. The recovery period will vary depending on the type of procedure, the type of implant, the patient's medical history and age, and other variables.

During this period, the area of abnormality surrounding the implant may also heal and recover. For example, if a cervical disc was implanted, the muscular structure surrounding this area, which may have previously been overcompensating because of the abnormality may now have returned to a normal state. These surrounding structure changes, may affect the way the implant was positioned within the patient or the performance char-

acteristics of the implant. In this regard, if the implant was positioned based upon abnormal surrounding structure, the implant may subsequently not provide the full range of motion as anticipated, thereby potentially resulting in further surgeries being required. Alternatively, the initially selected performance characteristics of the implant may have changed to due subsequent healing or other actions, thereby rendering the initial performance characteristics inappropriate for the current patient's condition. These performance characteristics can be any type of characteristics regarding the implant, including stiffness, actuation, loading, range of motion, etc. With the implant being an adjustable or tunable implant, corrections may be made to compensate for any subsequent anomalies observed by the surgeon. Again, the anomalies may result from healing of surrounding tissue, incorrect initial placement, changes in performance characteristics, or any other reasons. It should also be pointed out that if undesirable performance characteristics result after healing, the surrounding tissue and discs may also be damaged or deteriorate, thereby compounding recovery time and maybe requiring additional implants. This is the reason that providing the proper performance characteristics after healing is so critical.

[0105] After the patient has healed for some time, a post-operative exam is performed, exemplified at block 248. This post-operative exam may be conducted in different manners, depending upon the type of implant, the type of sensors and controls available with the implant, as well as the types of adjustments available with the implant. Some implants may have adjustment capabilities that require minimally invasive percutaneous type procedures, while other implants can be adjusted telemetrically or adaptively, as further discussed herein. The pre-operative exam may also be carried out using various types of equipment, again depending upon the capabilities of the implanted device, further discussed herein.

[0106] The pre-operative exam includes a motion analysis study, represented by block 250. This motion analysis study generally involves articulating the area of interest to determine range of motion, strength, etc. During this motion analysis study, the patient 14 is typically put through various motion testing. This testing may include various calisthenics, treadmill performance, weight lifting, gate analysis, etc. The motion analysis 250 can be performed and studied using an imagebased procedure, set out at block 252, a sensor-based procedure, set out at block 254, or an image and sensorbased procedure, set out in block 256. It should also be pointed out that while block 250 is labeled motion analysis, the analysis can be performed via static imagebased procedures or static sensor-based procedures, which are contemplated and included in the motion analysis study 250. In this regard, as opposed to putting the patient through various motion tests, the static image data or sensed data can be obtained and reviewed, via the image-based block 252 or the sensor-based block 254 to determine if the performance characteristics have changed. These static studies would simply look at the proper placement, impingement, etc. in the areas of interest to be used for subsequent post-operative tuning, further discussed herein.

[0107] The image-based procedure may be performed by either employing a localization or navigation tracking system or capturing image data, such as 3D or 4D image data, by an imaging device, such as a 4D CT imaging device. Should the motion analysis study be performed using localization or navigation technology, capturing image data and registration is performed as disclosed herein. U.S. Patent No. 6,434,415, entitled "System for Use in Displaying Images of a Body Part," issued August 13, 2002, also discloses pre-operative planning using navigation technology, which is hereby incorporated by reference. In general, pre-acquired image data may be obtained, for example, in the cervical spinal region. Before this image data is obtained, fiducial markers and localization sensors may be attached to each vertebrae of interest. Once the image data has been captured with these sensors in place, the patient 14 may be positioned on a treadmill with the tracking system 44 placed in proximity to track the motion of each vertebrae. This motion can include a gate analysis study of the patient's motion as well. Before the motion analysis begins, the navigation space of the patient 14 is registered to the pre-acquired images. Once the patient 14 begins the motion or movement for the motion analysis 250, tracking of the moving vertebrae may be captured and illustrated on a display, such as the display 10, or any other display.

[0108] If localization and navigation technology is not employed, image data may simply be captured over time during the motion analysis 250, for example, by the use of a four-dimensional CT scan. With this image data .. captured, each individual vertebrae may be segmented out using known segmenting algorithms. These types of algorithms generally involve thresholding or templates, which will segment out each vertebra in the scan. Once each vertebrae is segmented out, finite element analysis may be performed using known finite element analysis. The finite element analysis may also be used to calculate the instantaneous center of rotation "x". The information gathered during motion analysis 250 is used to determine the necessary adjustment of the implant at block 258. This information may include visualization of impinged areas around the implant, misalignment, etc. [0109] Should the motion analysis be sensor-based, as illustrated at block 254, the sensor readings of various parameters are used to determine if there is any necessary adjustment, at block 258. The sensor based approach may either take readings from sensors located within the implant or from sensors attached to the patient during this analysis. The sensors may take temperature readings, which can indicate potential friction and higher forces, strain or stress readings, as well as load readings or any other parameter readings. Again,

this information is used at block 258 to determine the necessary adjustment to the implant.

[0110] At block 256, both an image and sensor-based motion analysis may be conducted. This analysis essentially combines the image data at block 252 and the sensor data at block 254 to perform the post-operative analysis of the patient. Again, this information is used at block 258 to determine any necessary adjustments of the implant. When using both the image and sensor-based motion analysis, the sensed parameters may be synchronized in time with the image data to provide information on when the sensed parameters were captured relative to the time and the image.

[0111] At block 258, the data captured during motion analysis 250 is studied to determine whether any adjustments are necessary relative to the implant. For example, if a cervical disc was implanted and the patient healed and subsequent spinal alignment occurred, the range of motion may be compromised. In order to provide the proper range of motion, post-operative tuning of the implant may be necessary, based on the motion analysis study 250.

[0112] The post-operative tuning of the implant may also be necessary when the performance characteristics of the implant have changed.

Performance characteristics may be selected, based on various criteria, such as when the patient is in a relatively static state, thus requiring certain performance characteristics, as compared to when the patient is in vigorous active state, where the performance characteristics must be changed. For example, the spinal implant may not need significant stiffness in a relatively static condition, while in very active condition, the spinal implant may require a stiffer cushioning. The performance characteristic may have been selected when the patient was disabled, so that once the patient heals, the performance characteristics may have to be adjusted accordingly. This adjustment may be conducted using a minimally invasive adjustment procedure at block 260 or a telemetric adjustment procedure at block 262.

[0113] In the minimally invasive adjustment at block 260, percutaneous adjustment of the implant may be performed by actuating various adjustment mechanisms within the implant, further discussed herein. For example, adjustment screws may be positioned at hinge points within the implant and engaged by a driver in a minimally invasive type procedure to provide the proper adjustment, thereby reacquiring the proper range of motion, via adjusting the articulating surfaces of the implant. Adjustment of the performance characteristic, such as stiffness may also be performed, as further discussed herein.

[0114] Should a telemetric adjustment procedure be performed at block 262, a non-surgical adjustment would be performed. In this regard, the implant may be driven telemetrically, using known telemetric type wireless systems, such as that disclosed in U.S. Patent No. 6,474,341, entitled "Surgical Communication Power

System," issued November 5, 2002, which is hereby incorporated by reference or any other known wireless telemetric systems. The telemetric system may be an RF based or electromagnetic based telemetric system, as is known in the art. The implant may be a passive or active battery powered device that includes motors, pumps or any other devices used to adjust the implant, further discussed herein.

[0115] Once the adjustments have been performed, the procedure proceeds to block 264 where the adjustment is confirmed. If the adjustment is proper, the procedure ends at block 266. If not, further adjustments are performed. This pre-operative and post-operative procedure provides better initial implantation accuracy and implant selection, as well as the opportunity for post-operative tuning or adjustment of the implant. The post-operative tuning enables adjustment of articulating surfaces, supports, or other parameters within the implant post-operatively without requiring revision surgery.

[0116] Referring to Figures 11 and 12, an instrument assembly that includes a mounting platform 268 and an attachment jig 270 for use in a surgical navigated spinal procedure is illustrated. The mounting platform 268 is percutaneously attached to a series of vertebrae 272, via multiple K-wires 274.

In this regard, the mounting platform 268 is positioned outside the patient's body and above the vertebrae 272 of interest. The mounting platform 268 may be sized to span any number of vertebrae 272. In this example, four vertebrae 272 are spanned with the mounting platform 268 with the first and fourth vertebrae being secured to the mounting platform, via K-wires 274. The mounting platform 268 is designed to retain the spanned or captured vertebrae 272 in a relatively fixed or rigid manner during the spinal disc implant procedure. With the first and last vertebrae captured via the K-wires, the intermediate vertebrae 272 are generally held in a substantially fixed manner. Upon removing cartilage and other intermediate material between adjacent vertebrae, additional K-wires may be necessary for the intermediate vertebrae 272 to maintain the rigid structure.

[0117] The mounting platform 268 generally includes a rectangular-shaped beam 276 and a pair of outer attachment members 278. The rectangular beam 276 defines a plurality of peg holes 280, which are used to adjustably and removably retain the jig 270, along the member 276. The rectangular beam 276 also defines access and viewing holes or ports 282 enabling access from above and viewing of the relevant vertebrae. These access windows 282 can also be used to receive or pass surgical instruments during the medical procedure. Each attachment member 276 defines K-wire holes 284, which slidably receive the K-wires 274 in order to retain and secure the mounting platform 268 relative to the vertebrae 272.

[0118] An exemplary positioning jig 270 is illustrated in further detail in Figure 12 and is operable to be removably attached to the mounting platform 268, as

shown in Figure 11. In this regard, the jig 270 includes attachment peg 286 that is slidably received within holes 280. Positioned adjacent to the peg 286 is a pair of shoulders 288 that extend on either side of the rectangular beam 276 as the peg 286 is received within the hole 280. The jig 270 may be positioned along any part of the rectangular beam 276 by simply slidably inserting the peg 286 into one of the selected holes 280. Alternatively, any type of attachment mechanism to attach the jig 270 to the mounting platform 268 may be used. Once the peg 286 is positioned in one of the selected holes 280, the jig 270 is positioned substantially between a pair of vertebrae 272 in which the surgical procedure will be performed. The jig 270 further includes a work platform 290 that defines a passage 292 and includes a securing mechanism 294. The work platform 290 is positioned at an angle relative to the mounting platform 268 to provide intervertebral body access. The angle 296 illustrated with jig 270 provides for a 45° working platform 290. It should also be pointed out that multiple jigs 270 may be provided with the working platform 290 being positioned at various angles or the jig 270 may be adjustable to vary the angle, via a hinge or an adjustment mechanism between the working platform and the body of the jig 270.

[0119] The working platform 290 enables various instruments to be attached to the working platform, via the attachment mechanism 294, which may be a screw attachment, quick lock attachment, snap-fit attachment, or any other type of attachment mechanism. In one embodiment, a robot 298 may be attached to the working platform 290. This robot 298 may be remotely controlled and be used to drive milling, drilling, resection, or other instruments 300 through the passage 292. The robot 298 can either actuate the motor for the instrument 300 or can simply provide and act as an adjustable guide tube that may be controlled directly or remotely. Any type of known robotically controlled instrument may be utilized. Alternatively, the jig 270 may retain a manually adjustable guide tube that receives various instruments to be used during the procedure. The adjustable guide tube may also be lockable into a desired position in order to provide a rigid guide tube. Still further, the jig 270 may simply be used to pass and guide various instruments between the vertebral bodies 272. In this regard, the instruments as illustrated in Figures 8a-8g may be used in accordance with the jig 270 or other jigs providing various types of access ports 292. These access ports may be circular, slotted or any other shaped port to enable access between the vertebral bodies 272.

[0120] Generally, the mounting platform 268 will not include any localization sensors 58 or fiducial markers 60. The localization sensors 58 are generally positioned relative to the jig 270. The localization sensors 58 may be positioned on the guide tube and on the surgical instrument to determine orientation and depth of the surgical instrument 300, respectively. The localization sensor 58 may also be positioned on the robotically control-

led device 298 to determine both orientation and depth of the instrument 300. The mounting platform 268 may also include localization sensors 58 if desired, which may be used to provide further localization of the vertebrae 272. It should further be pointed out that the dynamic reference frame 54 may be attached or integrated into the mounting platform 268 in order to provide increased accuracy during the implant procedure. In this regard, since any motion of the mounting platform 268 would be identified, via an integrated dynamic reference frame 54, this motion is positioned substantially adjacent to the area of interest and the area being operated upon, providing increased registration and tracking of the instruments during the procedure.

[0121] By providing the mounting platform 268 that spans multiple vertebrae, multiple segment implantation may be performed in a minimally invasive and surgical navigated manner between the multiple vertebrae 272. For example, as illustrated, three separate cervical discs may be positioned between the four vertebrae 272 without requiring removal or replacement of multiple jigs as would typically be necessary. By providing a mounting platform 268 that can accommodate various size jigs and can be positioned between various vertebrae 272, a more precise and accurate implantation may be achieved in a more minimally invasive and efficient manner.

A ball and socket type cervical disc implant 302 [0122] is illustrated in Figure 13 that provides for percutaneous adjustment. The cervical disc implant 302 is based upon the Prestige Cervical Disc System provided by Medtronic Sofamor Danek of Minneapolis. Minnesota, but includes a tuning or adjustment capability. It should also be pointed out that while a cervical disc implant is disclosed herein, the present invention is not limited to merely cervical disc implants, but may include thoracic and lumbar spinal implants, as well as any other type of orthopedic implant that may require post-operative tuning. The cervical disc 302 comprises two articulating members that include a socket member 304 and a ball member 306. The socket member 304 includes a mounting flange 308 that defines generally two mounting holes 310 for receiving bone screws 312. The socket member 304 also defines the articulating socket 314 and is generally placed at an angle relative to the flange 308. Located at the junction between the flange 308 and the socket 314 is an adjustment or hinge region 316 defining an adjustment slot 318. Located within the adjustment slot 318 is an adjustment screw 320. Upon percutaneously engaging a head 321 of the adjustment screw 320, via any known driving instrument, the angle 324 between the flange and the socket 314 may be adjusted, via a wedge portion 322, in a minimally invasive manner. The head 321 may include a hex, a Phillips, a slotted, or any other type of engagable drive mechanisms that can be engaged by any type of instrument. Moreover, the adjustment screw 320 may be reversed so that the head 321 is located opposite, at 90°, or at any other orientation other than as illustrated to provide a different access point for the head 321.

[0123] The ball member 306 also includes a flange 326 defining screw holes 328 to receive bone screws 312. The ball member 306 also includes an articulating ball or spherical surface 330 that articulates with the socket 314. The flange 306 also includes adjustment or tuning portion 332 that defines a slot 334 for receiving another set screw 320 having head 322. Again, upon adjustment of the set screw 320, the angle 336 between the flange 326 and the ball 330 is adjusted in a minimally invasive manner, via percutaneous placement of a surgical driver that engages the head 321 of the adjustment screw 320.

[0124] By providing tuning or adjustment portions 316 and 332 relative to the ball 330 and socket 314, adjustment of the articulating ball 330 relative to the socket 314 may be made. Again; after a motion analysis 250 has been performed, a minimally invasive adjustment of the implant 302, such as the implant shown in Figure 13 may be performed by simply adjusting set screws 320. This adjustment may relieve impingement, increase range of motion, or provide other post-operative adjustments that would previously require a revision type surgical procedure.

[0125] Referring now to Figure 14, a modified embodiment of the cervical disc 302 is illustrated. In this regard, like reference numerals will be used to identify like structures as shown in Figure 13. The implant 302 provides for a telemetric type adjustment, as well as telemetric sensing capabilities. In this regard, the socket member 304 includes an actuator/controller 338 and a sensor 340 positioned along the articulating surface of the socket 314. Likewise, the ball member 306 also includes an actuator/controller 342 and a sensor 344 positioned along the articulating ball surface 330. The sensors 340 and 344 may be used to sense various parameters in the articulating joint, including temperature, pressure, stresses, strain and other loading properties. These sensors 340 and 344 may be used during the sensor based motion analysis 254 to sense the noted parameters during the motion analysis study 250. This sensed information is sent to its corresponding actuator/controller 338 or 342, which is able to telemetrically transmit information, further discussed herein, during this sensor based motion analysis 354.

[0126] Each actuator/controller 338 and 342 may either be a passive type device or an active rechargeable battery powered device. If the actuator/controllers 338 and 342 are passive type devices, they may include resonant LC circuits, which will resonate when adjacent generating coils, generate an electromagnetic field, thereby enabling transmission of the sensed information from sensors 340 and 344. An example of such a system is set out in U.S. Patent No. 6,474,341, entitled "Surgical Communication and Power System," issued November 5, 2002, which is hereby incorporated by reference. Other types of known wireless telemetric systems may also

be utilized. Actuator/controllers 338 and 342 may also be battery powered using rechargeable batteries that are either embedded within the implant or positioned remote from the implant and implanted within the patient, similar to known pacemaker technology. These rechargeable batteries may be recharged telemetrically similar to existing pacemaker batteries, as is known in the art.

[0127] If the system is a passive system, the data may be acquired from the corresponding sensor during the motion analysis study 250 in the post-operative exam 248 during the various motion tests performed on the patient 14. This information is gathered at the time of the study and is used to analyze whether or not further adjustments are necessary to the implant 302. Alternatively, if the system is an active system and battery powered, data may be sampled over time, stored in memory and transferred during the motion analysis study 250 or during other transfer periods, as further discussed herein. With this type of telemetric system, the implant 302 may be adjusted remotely by driving either actuator/controller 338 or 342 to remotely adjust the adjustable set screw 320, via known actuation type mechanisms. Again, while a hinge/set screw adjustment mechanism is shown, any other appropriate adjustment mechanism may be employed, such as worm gears, pinions, etc. Thus, telemetric adjustment 262 may be performed by simply positioning a corresponding transmit and receiving instrument adjacent to the implant site to both receive sensor information and remotely drive the actuators/controllers 338 and 342 to provide remote telemetric adjustment of the implant in a non-surgical manner. By adjusting either the angle 324 or the angle 336, the range of motion, contact, articulating surface adjustments, or other type of adjustments to relieve impingement and increase the range of motion may be performed in a post tuning technique. Briefly, Figure 15 shows the implant 302 implanted between a pair of vertebrae 272 and aligned, such that the instantaneous center of rotation are properly positioned within the center articulating longitudinal axis Y of the spine (see Fig. 8a).

[0128] Referring to Figure 16, another embodiment of a cervical spinal implant 346 is illustrated. The spinal implant 346 is based on the spinal disc prosthesis, set out in U.S. Patent No. 5,674,296, entitled "Human Spinal Disc Prosthesis," issued October 7, 1997 and U.S. Patent No. 5,865,846, entitled "Human Spinal Disc Prosthesis," issued February 2, 1999, each of which are hereby incorporated by reference. That is also known as the "Bryan Cervical Disc System," offered by Medtronic Sofamor Danek of Minneapolis, Minnesota. The spinal implant 346, however, also includes a tuning and adjustment mechanism. The spinal implant includes a pair of rigid support plates 348 and a pair of attachment flanges 350 that define attachment holes to receive bone screws (see Figure 8g). Positioned between the support plates 348 is a flexible bladder device

352.

In order to provide for either minimally invasive [0129] or telemetric adjustment of the implant 346, the bladder mechanism 352 is separated into a plurality of individual bladders 354. As illustrated, the implant 346 includes three adjacent bladders 354. Located within each bladder 354 is a sensor 356 that is used to sense the pressure within each bladder 354. These sensor readings are passed to a bladder control system 358. The bladder control system 358 may again be a passive device or an active battery powered device. If passive, the sensor information will be received during the motion analysis study 250 and adjustment may be performed telemetrically during this study using known telemetric driving devices. If the bladder control system 358 is an active powered system, the system may either operate similar to the passive system or may be an adaptive system that provides real time adjustment for the implant 346. In this regard, each sensor 356 may sense pressure differences in each bladder 354 while the bladder control system 358 attempts to equalize the pressures in the bladders 354 in a real time manner. The bladder control system 358 includes a processor controller 358 and either a battery or known passive driving device. The bladder control system 358 also includes a pump used to transfer fluid retained within the bladders 354 by controlling remote valves 360 and a memory if necessary for storing sampled data.

[0130] The implant 346 may also include a reservoir 361 that retains a drug that may be delivered through the external valve 360 and controlled by the bladder control system 358. In this way, controlled drug delivery to the surrounding bone may also be achieved with the implant 346. The drug can include a bone morphagenic protein (BMP) that is able to increase bone density and fusion of broken bones, by delivering the BMP over time to the surrounding infected bones. This drug delivery capability of the implant 346 may be actively delivered if the system is battery-powered, or telemetrically delivered, via an active or passive device during patient exams.

In operation, the implant 346 may be used to [0131] sense pressure in each individual bladder 354, via the sensors 356 during the post-operative motion analysis 250. With this information, a surgeon can direct the bladder control system 358 to compensate for any abnormalities in pressure in the bladders 354 in order to try to achieve uniform pressure throughout the implant 346. The bladders 354 generally include a saline solution that can be transferred between bladders 354, via the bladder control system 358 and control valves 360. In addition, there is an external valve 360 that may be used to release saline fluid harmlessly into the body to relieve pressure. Alternatively, the external valve 360 may be used to receive additional fluid percutaneously in a minimally invasive way. Thus, the implant 346 may be postoperatively adjusted or tuned, depending upon the healing of the patient, post-operative trauma, or to provide further refinement and increased performance of the implant 346.

[0132] Alternate embodiments of the implant 346 is illustrated in 16a-16c. Here again, like reference numerals are used to identify like structures. The spinal implant 346 is substantially similar to the spinal implant illustrated in Figure 16, except that the spinal implants illustrated in Figure 16a-16c are multi-segment implants 346 that allow for a minimally invasive technique and a posterior implantation approach. The implant 346 illustrated in Figure 16a includes a pair of rigid support plates 348 that include a hinged region 349. This hinged region 369 includes a single hinge that enables the implant 346 to be substantially compressed so that the plates 348 are adjacent to one another. Once adjacent to one another, the plates 348 may be folded via the hinge region 349 creating a semi-circular shape that is significantly smaller than the whole implant 346. This enables the implant to be implanted posteriorly in a minimally invasive manner by simply sliding the folded implant 346 into a small incision and re-assembling or unfolding the implant 346, along the hinge region 349 at the implant area. The hinge 349 also includes a lock 351 that is used to lock the hinge 349 to insure that each plate 348 is locked in a planar fashion. Once locked, the implant 346 is positioned between the adjacent vertebrae 242 similar to that shown in Figure 16.

[0133] Implant 346, illustrated in Figures 16b and 16c also includes a hinged region 349 that consists of a pair of hinges positioned on either side of the flange 350. Again, the hinge region 349 enables the end plates 348 to be folded, as illustrated in Figure 16c to enable a posterior minimally invasive procedure. This implant 346 also includes a lock 351 that rotates to lock the pair of hinges in the hinge region 349 in a substantially planar manner.

[0134] Another embodiment of the spinal implant 346 is shown in Figure 17, which provides a different type of adjustment mechanism. Here again, like reference numerals will be used to identify like structures. Again, the spinal implant 346 includes a pair of supporting plates 348, a pair of flanges 350 and a support or bladder device 352. Located within the bladder device 352 is a single bladder 354, which can be filled with a saline solution, or optionally not filled with fluid. Again, sensors 356 are located in different regions within the bladder 354 and used to either sense fluid pressure or used as a strain gauge to measure loading forces. The readings from the sensors 356 are read by a force control system 362, which can again either be a passive device or a battery powered active device. The force control system 362 operate similar to the bladder control system 358, except that as opposed to directing fluid between various bladder chambers, it includes force control beams or members 364 that are used to apply a force to the plurality of springs 366 positioned within the bladder 354. By compressing the springs 366 in different quadrants with the control beams 364, tension in the springs 366 are increased, thereby providing additional support within the implant 346. Each spring may be selectively adjusted, depending upon the desired tuning or adjustment necessary. Again, this adjustment is based upon the motion analysis study done during the post-operative exam 248.

[0135] The force control system 362 may be used to adaptively or actively adjust the implant 346 if the force control system is an active battery powered system. Alternatively, the force control system 362 may adjust the force within the implant 346 during the telemetric adjustment 262 if the system is simply passive. The bladder control system 358 and the force control system 362 may be formed using conventional micro electronics and mechanical devices or may be formed from micro electromechanical system (MEMS) technology, known in the art.

[0136] A multiple segment implantation is illustrated in Figure 18 that includes multiple implants 368. Each implant 368 may be implanted utilizing the mounting 🖺 platform 268 and jig 270, as illustrated in Figures 11 and 12. Implants 368 may also be implanted using other procedures, such as that shown in Figures 8a-8g. Each implant 368 includes a sensor 370 and an adjustment actuator 372, similar to that shown in Figure 14, 16, and 17. However, each implant 368 is controlled and actuated, via an active rechargeable battery powered external controller 374. Optionally, each implant 368 may include its own individual internal controller 374 that can communicate to the other implants 368, via a wireless or wire connection. Alternatively, a single internal master controller 374 may be positioned within one of the implants 368, which is used to control and drive the remaining implants 368 in a master/slave relationship.

[0137] Controller 374 is used to sense various parameters again, such as temperature, pressure, etc. where actuators 372 are used to tune or adjust each implant 368 accordingly. The controller 374 may be implanted adjacent to the spinal region, similar to a controller and battery for a pacemaker. The multiple segment implantation with each implant 368 communicating with the other surrounding implants 368 enable real time adaptive control of this spinal region, such as the cervical spi--- nal region of the patient 14. In other words, the controller 374 may sense, via the sensors 370 whether any one of the implants 368 is under too much pressure or one may be too laxed and adjust accordingly, depending upon the patient's movements. In this regard, when the patients at rest, extra support between the vertebrae 272 may not be necessary. However, when the patient 14 is doing physical activities or exercise, additional support may be necessary between each vertebrae 272 and each implant 368 may be expanded during this period in an adaptive manner. Alternatively, the controller may again simply be a passive controller or an active controller and used to send and receive information, as well as adjust the implants 368 during the post-operative exam 248, via the telemetric adjustment 262.

[0138] Turning to Figure 19, an exemplary telemetric system used for performing the motion analysis 250 is illustrated. In this regard, the patient 14 may undergo the motion analysis 250 by exercising on a treadmill 376. The treadmill 376 is positioned within a transmit/receive module 378. When the patient 14 is positioned within the transmit/receive module 378 and exercising on the treadmill 376, information can be collected from the particular implant during the motion analysis 250 using the sensor based 252 data analysis, via the telemetric adjustment 262. In other words, the transmitter/receive module 378 includes signal transmitters and receivers to either actuate a passive or active controller to receive sensed information. This information is forwarded to a control processor 380 where the surgeon can analyze... the collected sensed data. Once the data has been analyzed, the controller 380 is used to actuate the transmit/ receive module 378 to adjust one or more implants in the patient, via the control actuator circuits, disclosed above. It should also be noted that an imaging device may also be positioned adjacent to the patient 14 while the patient is on the treadmill 376 to provide both an image-based and a sensed-based motion analysis 250, as previously discussed.

[0139] Referring now to Figure 20, another telemetric system used to transmit motion analysis information to the doctor is disclosed. With this technique, the patient 14 can simply conduct a self analysis by positioning him or herself adjacent to a computer 382. Attached to the computer 382 is a transmit/receive module 384. The transmit receive module 384 operates similar to the transmit receive module 378, except that the patient 14 can simply run through a set of suggested motions, while the transmit receive module 382 telemetrically receives information from the implant positioned within the patient 14. This information can be transmitted, via the computer 382 online to a receiving hospital or doctor's office. The doctor may then analyze this information, make a recommendation to the patient 14 whether the patient 14 should come in to the office for a telemetric adjustment 362 of the patient's implant. Alternatively, the doctor may also simply instruct the transmit control module 384, via the computer 382, to perform the telemetric adjustment of the patient 14 in the patient's home.

[0140] The procedure 230, as well as the associated implants, systems and instruments, enables both preoperative and post-operative review and analysis. Additionally, post-operative tuning of the implant may also be achieved without requiring revision surgery or highly invasive types of procedures. In this regard, either minimally invasive or telemetric adjustments of the implants may be achieved.

[0141] The description of the invention is merely exemplary in nature and, thus, variations that do not depart from the gist of the invention are intended to be within the scope of the invention. Such variations are not to be regarded as a departure from the spirit and scope of the invention.

15

20

[0142] The method for tuning an implant positioned within a patient, can include conducting a post-operative exam of the implant positioned within the patient.

. ;

[0143] The post-operative exam can include a motion analysis on the operation of the implant positioned within the patient.

[0144] The motion analysis can be an image-based motion analysis and includes capturing images of the implant positioned within the patient.

[0145] The images can be captured during motion of the patient to observe the implant in operation.

[0146] The images can be captured during a static condition of the patient to observe the implant in operation.

[0147] The motion analysis can be a sensor-based study where various parameters are sensed at the implant site.

[0148] The parameters can be selected from a group comprising temperature, pressure, stress, strain, loading force, tortional force, and frictional force.

[0149] The sensed parameters can be sensed from at least one sensor positioned within the implant.

[0150] The motion analysis can be both image-based and sensor-based where the sensor-based analysis includes obtaining at least one reading of a parameter at the implant site.

[0151] The step of determining if any adjustment of the implant is necessary can include analyzing the motion analysis study.

[0152] The step of adjusting the implant can include performing a minimally invasive adjustment.

[0153] The minimally invasive adjustment can be performed by engaging an adjustment mechanism within the implant with an adjustment instrument.

[0154] The adjustment instrument can actuate an adjustment portion within the implant to adjust the orientation of the implant within the patient.

[0155] The adjustment instrument can engage the implant to at least deliver a fluid or release a fluid within the implant.

[0156] The step of adjusting the implant can includes adjusting the implant telemetrically in a non-surgical procedure.

[0157] The telemetric adjustment can include wirelessly actuating a driver within the implant to adjust performance of the implant positioned within the patient.

[0158] The wireless actuation of the driver further can include wirelessly activating at least one of a passive driver and an active driver positioned within the implant.

[0159] The step of adjusting the implant can be a realtime adaptive adjustment of the implant without any outside intervention.

[0160] The method can further comprise tuning a plurality of implants positioned within the patient.

[0161] The method can further comprise adaptively adjusting the plurality of implants with a single controller positioned within the patient.

[0162] The method can further comprise positioning

the controller within one of the plurality of implants.

[0163] The method can further comprise providing communication between the plurality of implants.

[0164] The implant used in the method can be a spinal disc implant, e.g. a cervical spinal disc implant, and can include an articulating ball and socket and/or an adjustable bladder.

[0165] A further aspect of the invention relates to a display for use in guiding a medical device to a target in a patient during a medical procedure, said display comprising:

a first indicia illustrated on said display identifying a first degree of freedom information; a second indicia illustrated on said display identifying a second degree of freedom information; a third indicia illustrated on said display identifying a third degree of freedom information; a fourth indicia illustrated on said display identifying a fourth degree of freedom information; and a fifth indicia illustrated on said display identifying a fifth degree of freedom information, wherein said first, second, third, fourth, and fifth indicias are illustrated on said display to assist in guiding the med-

[0166] The display can further comprising a sixth indicia illustrated on said display identifying a sixth degree of freedom of information.

ical device to the target.

[0167] The indicia can be illustrated in a coordinate system selected from a group comprising a Cartesian coordinate system, a spherical coordinate system, a polar coordinate system, and a combination thereof.

[0168] The first, second, and third indicia can be X, Y, and Z axes, respectively, and the fourth, fifth, and sixth indicia can be rotation about the X, Y, and Z axes, respectively.

[0169] At least one indicia can illustrate a safe zone in relation to said corresponding degree of freedom information.

[0170] Each indicia can illustrate a safe zone.

[0171] The safe zone can be illustrated by way of color, audible tones, cross-hatching and a combination thereof.

[0172] The display can further comprise an indicia of a tip portion of the medical device and an indicia of a hind portion of the medical device illustrated on said display depicting the relative position and orientation of the medical device.

[0173] This tip can be represented by cross-hairs and the hind is represented by a circle.

[0174] Said X axis can be represented by a horizontal axis, said Y axis by a vertical axis crossing said X axis, and said Z axis by a separate axis having an arrow identifying the position along the Z axis.

[0175] Rotation about the X axis can be illustrated by a straight bar graph having an arrow identifying the amount of rotation, wherein the rotation about the Y axis

is illustrated by a straight bar graph having an arrow illustrating the amount of rotation about the Y axis and the rotation about the Z axis is illustrated by an arcuate bar graph having an arrow illustrating the rotation about the Z axis.

[0176] The medical procedure can be selected from an orthopedic procedure, a spinal procedure, a neurovascular procedure, a cardiovascular procedure, a soft tissue procedure, and a combination thereof.

[0177] The medical device can be selected from a group comprising a medical instrument, a surgical implant, and a combination thereof.

[0178] The medical instrument can be selected from a group comprising a catheter, a reamer, a drill, a cutting block, a saw, an impactor, a biopsy needle, a delivery needle, and a combination thereof.

[0179] The implant can be selected from a group comprising an orthopedic implant, a spinal implant, a cardiovascular implant, a neurovascular implant, a soft tissue implant, and a combination thereof.

[0180] The display can further comprise a target/trajectory indicia identifying the target/trajectory for the medical device.

[0181] The display can further comprise indicia of a second target/trajectory identifying a second target/trajectory.

[0182] The display can further comprise at least one of an image of the patient superimposed on said display and one of a model of the patient superimposed on said display.

[0183] The display can further include multiple displays.

[0184] The first, second, third, fourth and fifth indicia can be simultaneously illustrated on said display.

[0185] The display can assist in the control of a surgical robot used during the medical procedure.

[0186] The target indicia can be an origin.

[0187] A further aspect of the invention relates to a navigation system for use in guiding a medical device to a target in a patient during a medical procedure, said navigation system comprising:

a tracking sensor associated with the medical device and operable to be used to track the medical device;

a tracking device operable to track the medical device with the tracking sensor; and

a display having indicia illustrating at least five degrees of freedom information and indicia of the medical device in relation to said at least five degree of freedom information.

[0188] The tracking sensor can be attached to the medical device.

[0189] The tracking sensor can be removably attached to the medical device.

[0190] The tracking sensor can be an electromagnetic tracking sensor.

[0191] The tracking sensor can be selected from a group comprising an electromagnetic sensor, an optical sensor, an optical reflector, an electromagnetic transmitter, a conductive sensor, an acoustic sensor, and combinations thereof.

[0192] The tracking device can be an electromagnetic tracking device.

[0193] The tracking device can be selected from a group comprising an electromagnetic tracking device, an optical tracking device, a conductive tracking device, an acoustic tracking device, and combinations thereof.

[0194] The navigation system can further comprise an imaging device operable to generate image data of a region of the patient.

[0195] The imaging device can be selected from a group comprising a C-arm fluoroscopic imager, a magnetic resonance imager (MRI), a computed tomography (CT) imager, and a positron emission tomography (PET) imager, an isocentric fluoroscopy imager, a bi-plane fluoroscopy imager, an ultrasound imager, a multi-slice computed tomography (MSCT) imager, a high-frequency ultrasound (HIFU) imager, an optical coherence tomography (OCT) imager, an intra-vascular ultrasound imager (IVUS), a 2D, 3D or 4D ultrasound imager, an intra-operative CT imager, an intra-operative MRI imager, and a single photon emission computer tomography (SPECT) imager.

[0196] The navigation system can be used as an image-less based navigation system wherein said image-less navigation system is operable to generate a model representing a region of the patient.

[0197] The navigation system can further comprise a dynamic reference frame attached to the patient and operable to be used to track movement of the patient relative to said tracking device.

[0198] The medical device can be a medical instrument or implant.

[0199] The medical instrument can be selected from a group comprising a catheter, a reamer, a drill, a cutting block, a saw, an impactor, a biopsy needle, a delivery needle, and combinations thereof.

[0200] The implant can be selected from a group comprising an orthopedic implant, a spinal implant, a cardiovascular implant, a neurovascular implant, a soft tissue implant, or a combination thereof.

[0201] The display can further include indicia illustrating a sixth degree of freedom information.

[0202] The first, second, and third indicia can be X, Y, and Z axes, respectively, and the fourth, fifth, and sixth indicia can be rotation about the X, Y, and Z axes, respectively.

[0203] At least one indicia can illustrates a safe zone in relation to said corresponding degree of freedom information.

[0204] The indicia of the medical device can include indicia identifying a tip portion and indicia identifying a hind portion of the medical device.

[0205] The tip can include indicia of crosshairs and

said hind includes indicia of a circle.

[0206] The medical device can have at least a first component and a second component where placement of said first component is dependent upon placement of said second component and wherein said display is operable to display proper position of said first component relative to said second component.

[0207] The medical procedure can be selected from an orthopedic procedure, a spinal procedure, a neurovascular procedure, a cardiovascular procedure, a soft tissue procedure, and a combination thereof.

[0208] A further aspect of the invention relates to a method for navigating and displaying a medical device during a medical procedure, said method comprising:

selecting a target for navigating the medical device to;

displaying the target in a coordinate system; tracking the medical device; and displaying the medical device in relation to the target with at least five degrees of freedom information.

[0209] The method can further comprise displaying indicia of the at least five degree of freedom information. [0210] The method can further comprise selecting a safe zone for at least one of the degrees of freedom coordinate.

[0211] The method can further comprise selecting a safe zone for each degree of freedom coordinate.

[0212] The method can further comprise selecting a type of coordinate system.

[0213] The method can further comprise aligning the medical device with the target using the five degree of freedom information.

[0214] The method can further comprise selecting a second target that is positioned in relation to the first target for a second medical device; and displaying the second target in the coordinate system.

The method can further comprise: [0215]

tracking the second medical device; and displaying the second medical device in relation to the second target with the at least five degree of freedom information.

The method can further comprise aligning the [0216] second medical instrument with the second target.

[0217] The method can further comprise:

determining if the alignment of the first medical device and the second medical device are correct; and realigning the first and/or medical device if necessary.

[0218] The step of selecting the target can further include selecting a center of a joint.

[0219] The step of selecting the target can further include selecting a target based upon an articulating plane of a joint.

[0220] The step of selecting the target can further include selecting an area between a first vertebra and a second vertebra for placement of an implant.

[0221] The step of selecting the target can further include selecting a site within an artery.

# Claims

15

A tunable implant positioned within a patient, said tunable implant comprising:

> a securing mechanism operable to secure said implant within the patient; an actuation portion operable to permit said tunable implant to move; and an adjustable portion operable to permit adjust-

> ment of said tunable implant after said tunable implant is positioned within the patient, wherein upon adjusting the implant, the implant provides improved performance.

- The tunable implant as defined in Claim 1 wherein said securing mechanism includes a flange operable to engage a bone of the patient.
- The tunable implant as defined in Claim 2 wherein *30* 3. said flange defines a passageway operable to receive a bone screw to secure the tunable implant to the bone.
- The tunable implant as defined in Claim 2 or 3, fur-35 ther comprising a pair of flanges operable to secure a first portion and a second portion of said tunable implant to a pair of articulating bones.
- The tunable implant as defined in any one of the preceding claims, wherein said actuation portion includes a first articulating surface and a second articulating surface, wherein said first articulating surface is operable to articulate with said second articulating surface. 45
  - 6. The tunable implant as defined in Claim 5, wherein said first articulating surface is a ball and said second articulating surface is a socket, wherein said ball articulates within said socket.
  - 7. The tunable implant as defined in any one of the preceding claims, wherein said actuation portion includes a bladder portion operable to permit said tunable implant to move relative to bone.
  - 8. The tunable implant as defined in Claim 7, wherein said bladder portion is filled with fluid.

24

50

*55* 

BNSDOCID. < EP\_\_\_\_\_1442715A2\_I\_>

10

15

20

25

30

35

40

45

- 9. The tunable implant as defined in Claim 7 or 8, wherein said bladder portion includes a plurality of bladders.
- 10. The tunable implant as defined in Claim 9, wherein said actuation portion includes a bladder control system operable to adjust pressure within each of said plurality of bladders.
- 11. The tunable implant as defined in Claim 10, wherein said bladder control system adaptively controls pressure in each of said plurality of bladders.
- 12. The tunable implant as defined in one of Claims 7 9, wherein the bladder includes at least one mechanical adjustment mechanism.
- 13. The tunable implant as defined in Claim 12, further comprising a force control system operable to control the mechanical adjustment mechanism within the said bladder.
- 14. The tunable implant as defined in Claim 10 or 13, further comprising at least one sensor operable to sense a parameter within said/each bladder.
- 15. The tunable implant as defined in Claim 14, wherein said control system is further operable to receive sensor readings from each of said sensors.
- 16. The tunable implant as defined in claim 10 or 13. wherein said control system is one of at least a passive and an active control system.
- 17. The tunable implant as defined in Claim 10 or 13, wherein said control system is operable to be controlled telemetrically in a non-surgical procedure.
- 18. The tunable implant as defined in any one of the preceding claims, further comprising at least one sensor positioned within said implant and operable to sense a parameter.
- 19. The tunable implant as defined in Claim 18, wherein said parameter is selected from a group comprising temperature, pressure, stress, strain, loading force, tortional force, and frictional force.
- 20. The tunable implant as defined in Claim 18 or 19, further comprising a controller operable to receive sensed information from said at least one sensor, said controller positioned within said implant.
- 21. The tunable implant as defined in any one of the preceding claims, wherein said actuation portion is a mechanical actuation portion operable to be actuated by an instrument during a minimally invasive procedure.

- 22. The tunable implant as defined in Claim 21, wherein said actuation portion includes a hinged region and an adjustable screw operable to be engaged by said instrument during the minimally invasive procedure.
- 23. The tunable implant as defined in any one of the preceding claims, wherein said actuation portion is an actuator operable to be actuated telemetrically in a non-surgical procedure.
- 24. The tunable implant as defined in Claim 23, wherein said actuation portion includes a hinged region and an adjustable screw wherein said actuator is operable to telemetrically adjust said adjustable screw during the non-surgical procedure.
- 25. The tunable implant as defined in any one of the preceding claims, further comprising a plurality of tunable implants, each including said securing mechanism, said actuation portion and said adjustment portion.
- 26. The tunable implant as defined in Claim 25, further comprising at least one controller operable to control adjustment in each of said plurality of tunable implants.
- 27. The tunable implant as defined in Claim 26, further comprising at least one sensor positioned in each implant operable to forward sensed information to said at least one controller.
- 28. The tunable implant as defined in Claim 26 or 27, wherein said at least one controller is operable to adaptively control each of said plurality of implants.
- 29. The tunable implant as defined in Claim 26 28, wherein said at least one controller is located within at least one implant.
- 30. The tunable implant as defined in Claim 26 28, wherein said at least one controller is remote from each of said implants and positioned within the patient.
- 31. The tunable implant as defined in Claim 26 30, wherein said at least one controller is operable to be controlled telemetrically.
- 32. The tunable implant as defined in any one of the preceding claims, further comprising a segmented portion operable to enable said implant to be folded in a manner to make the implant substantially smaller than its unfolded state.
  - 33. A tunable implant system for use in adjusting a tunable implant in a patient, said tunable implant system comprising:

20

30

*35* 

40

50

a tunable implant having an adjustment portion operable to permit adjustment of said tunable implant after said tunable implant is positioned within the patient; and a telemetric system operable to telemetrically receive data from said tunable implant, wherein said data is used to determine said adjustment

34. The tunable implant system as defined in Claim 33 wherein said tunable implant further includes an actuation portion operable to permit said tunable implant to move.

of said tunable implant.

- 35. The tunable implant system as defined in Claim 33 or 34 wherein said tunable implant further includes at least one sensor operable to sense a parameter within said implant, said parameter being at least part of said data received by said telemetric system.
- 36. The tunable implant system as defined in Claim 35 wherein said tunable implant further includes a controller operable to receive said sensed parameter information from said sensor and further operable to transmit said sensed information to said telemetric system.
- 37. The tunable implant system as defined in Claim 36 wherein said controller is at least one of a passive and an active controller.
- 38. The tunable implant system as defined in Claim 36 or 37 wherein said controller is further operable to actuate said adjustment portion to adjust said tunable implant.
- 39. The tunable implant system as defined in Claim 36, 37 or 38 wherein said telemetric system is operable to transmit and receive data from said controller.
- 40. The tunable implant system as defined in Claim 39 wherein said telemetric system is further operable to charge a battery in said controller.
- 41. The tunable implant system as defined in one of Claims 33 40 wherein said telemetric system further includes a controller operable to process said data received telemetrically by said telemetric system.
- 42. The telemetric implant system as defined in Claim 41 wherein said controller is further operable to transmit said data to a remote site.
- 43. An instrument assembly for use in performing spinal surgery, said instrument assembly comprising:

a mounting platform operable to be positioned

adjacent to a plurality of vertebrae; and a jig operable to be removably attached to said mounting platform, said jig further operable to support an instrument used during the spinal surgery.

- 44. The instrument assembly as defined in Claim 43 wherein said mounting platform is further operable to substantially rigidly secure said plurality of vertebrae.
- 45. The instrument assembly as defined in Claim 43 or 44 wherein said jig is further operable to be selectively and removably attached to said mounting platform, along a plurality of positions defined on said mounting platform.
- 46. The instrument assembly as defined in one of claims 43 45, wherein said instrument is a robotically controlled instrument coupled to said jig.
- 47. The instrument assembly as defined in one of claims 43 46, wherein a location of said instrument is tracked utilizing a surgical navigation system.
- 48. The instrument assembly as defined in Claim 47 further comprising a tracking device attached to said instrument assembly and operable to be tracked by said surgical navigation system.
- 49. The instrument assembly as defined in Claim 48 wherein said tracking device is attached to said jig.
- 50. The instrument assembly as defined in Claim 48 wherein said tracking device is attached to said instrument attached to said jig.
- 51. The instrument assembly as defined in one of claims 43 50, further comprising a plurality of jigs operable to be removably attached to said mounting platform, each of said jigs providing a different orientation of a working platform.
- 52. The instrument assembly as defined in Claim 51 wherein said jig further includes a working platform defining an opening for passing said instrument.
- 53. A method for tuning an implant positioned within a patient, said method comprising:

analyzing an operation of the implant positioned within the patient; determining if any adjustment of the implant is necessary; and adjusting the implant while the implant is positioned within the patient, wherein upon adjusting the implant, the implant provides improved performance.

**26** 

BNSUCCIU: <EP\_\_\_\_\_1442715A2\_1\_>





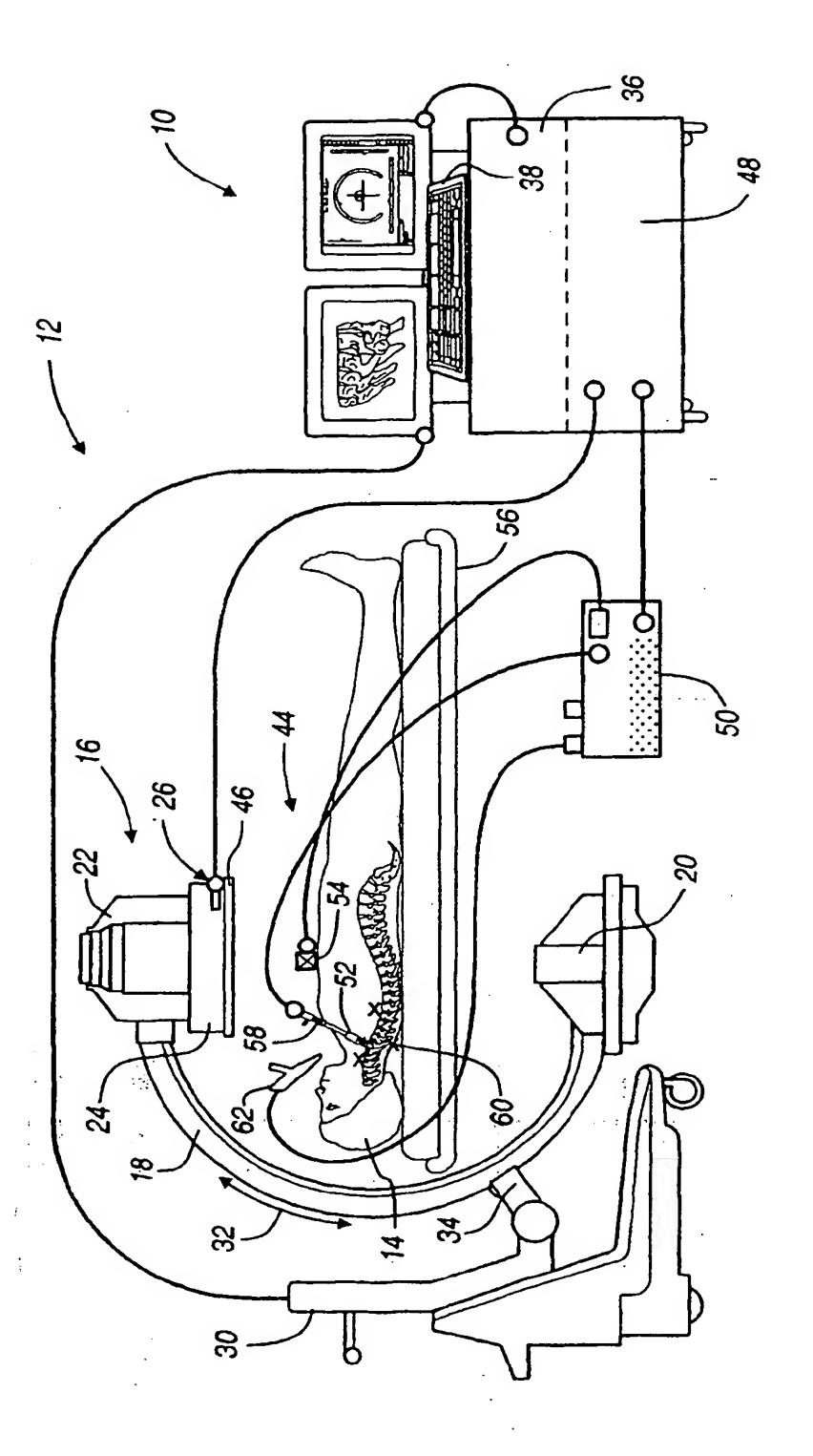



FIGURE 1

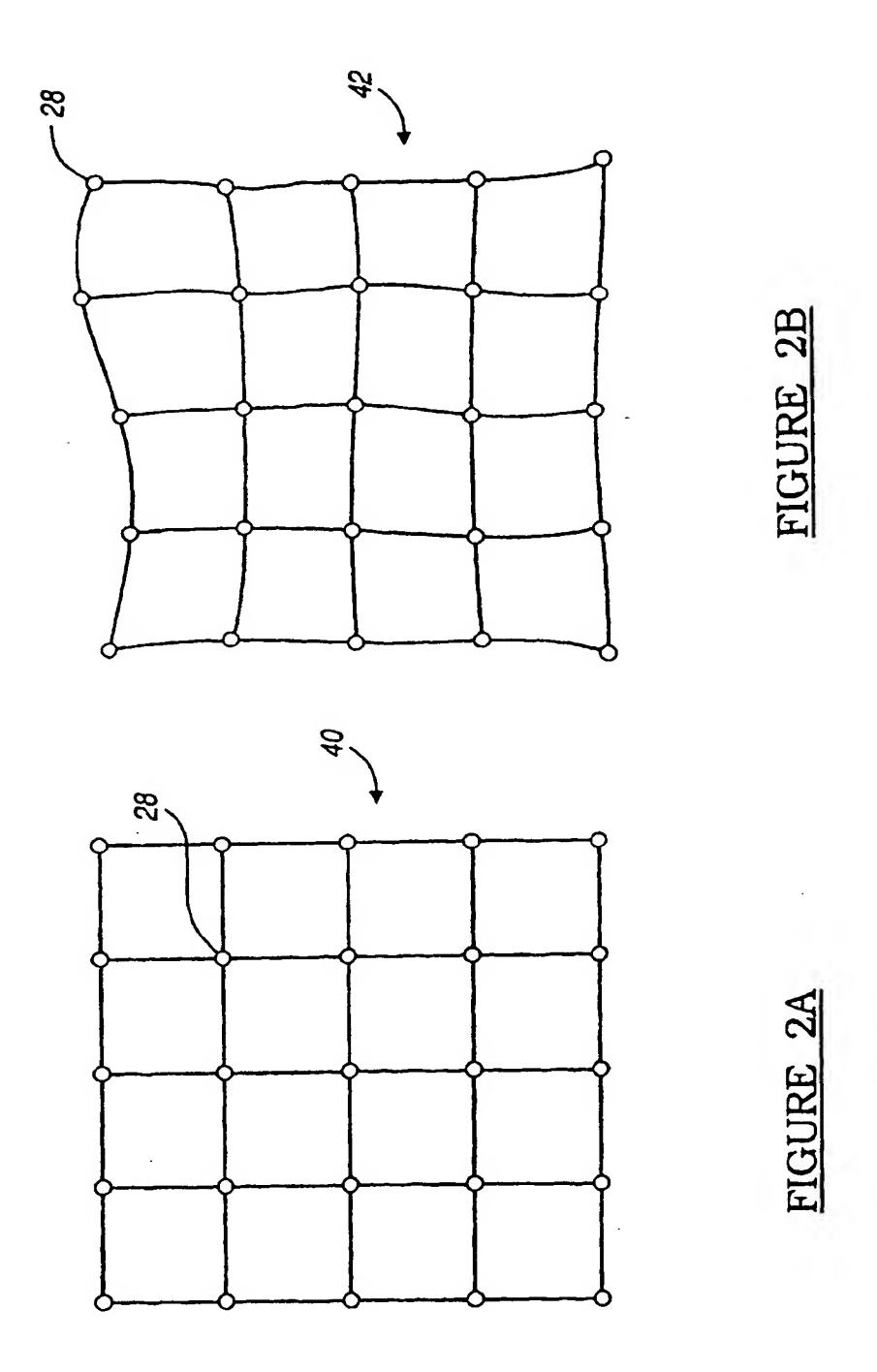


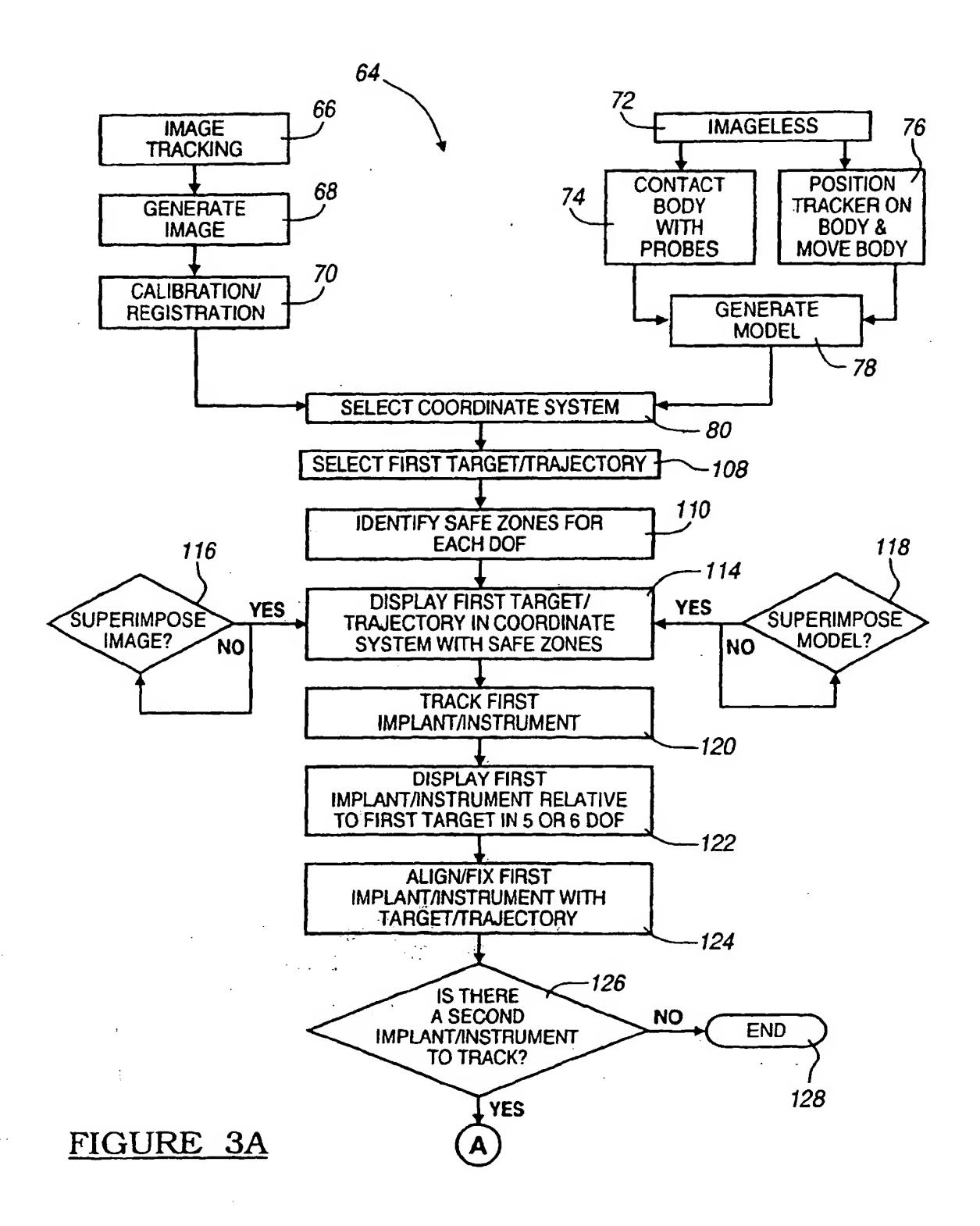


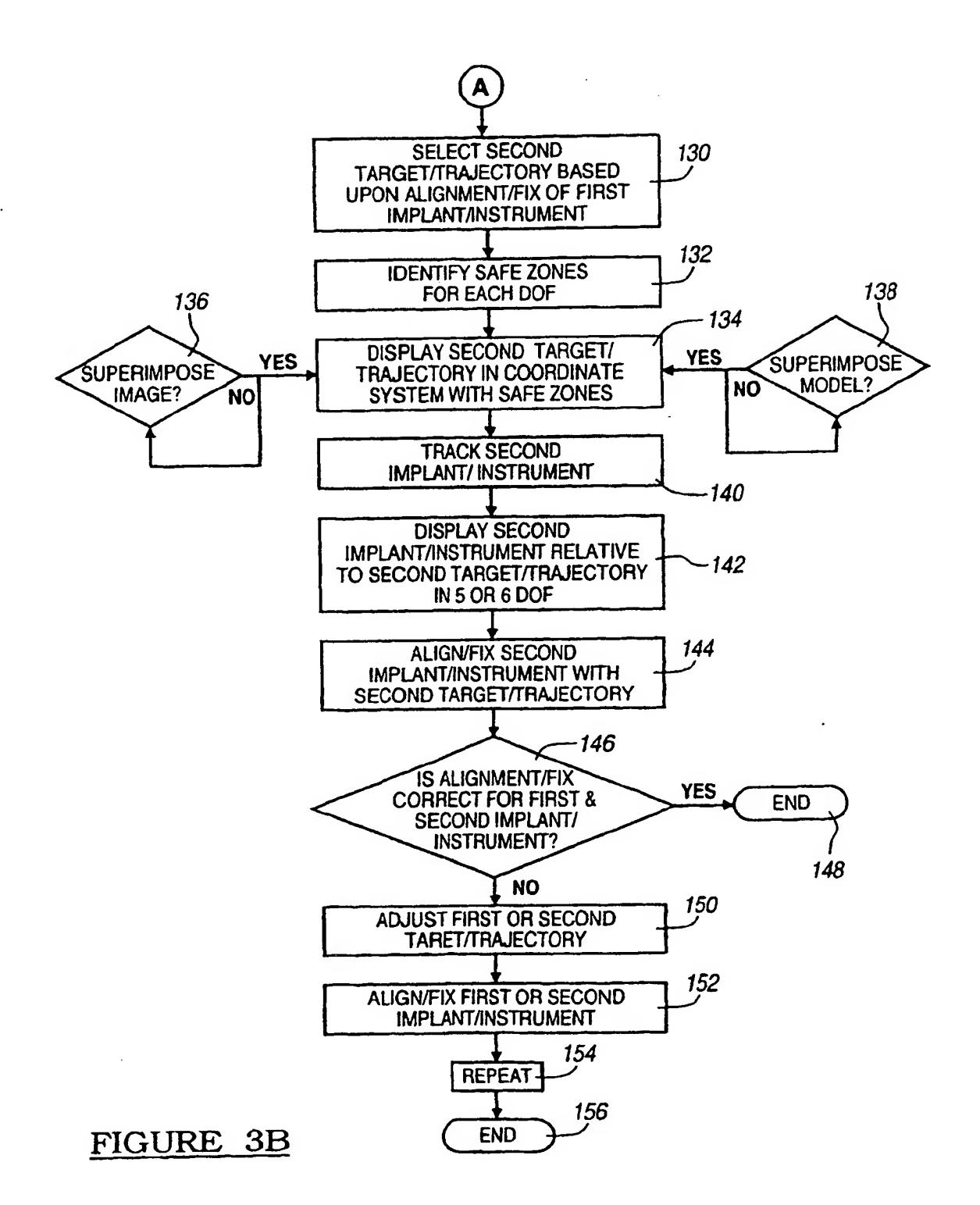


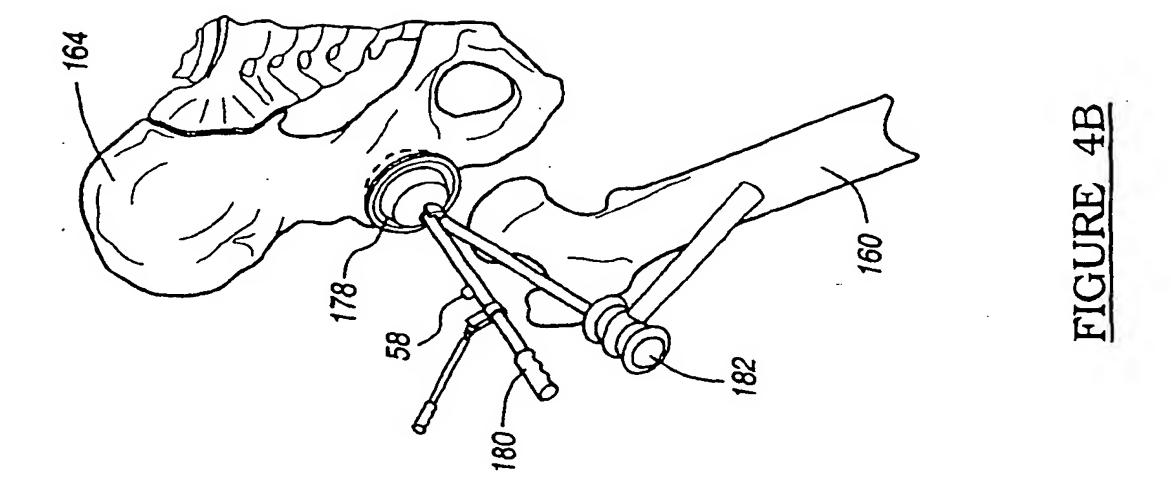


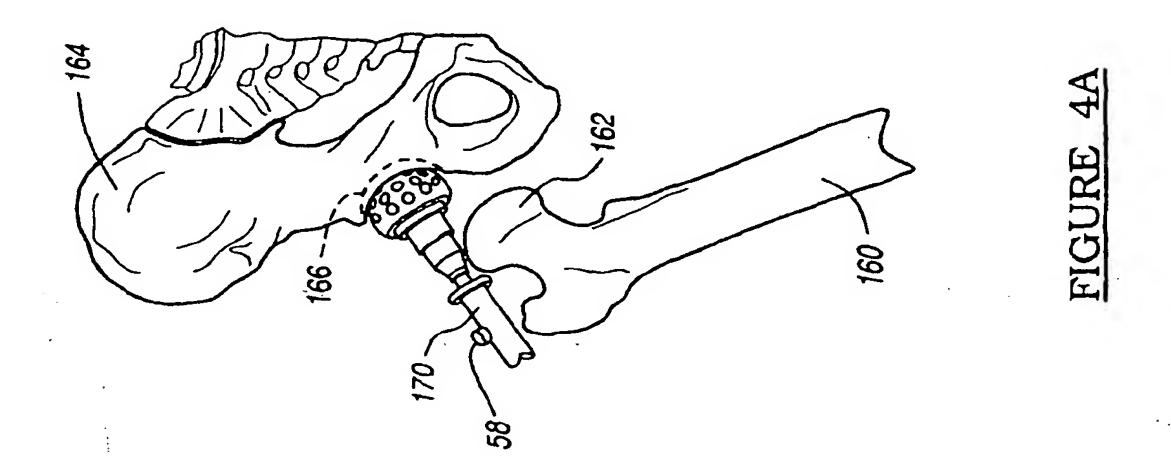

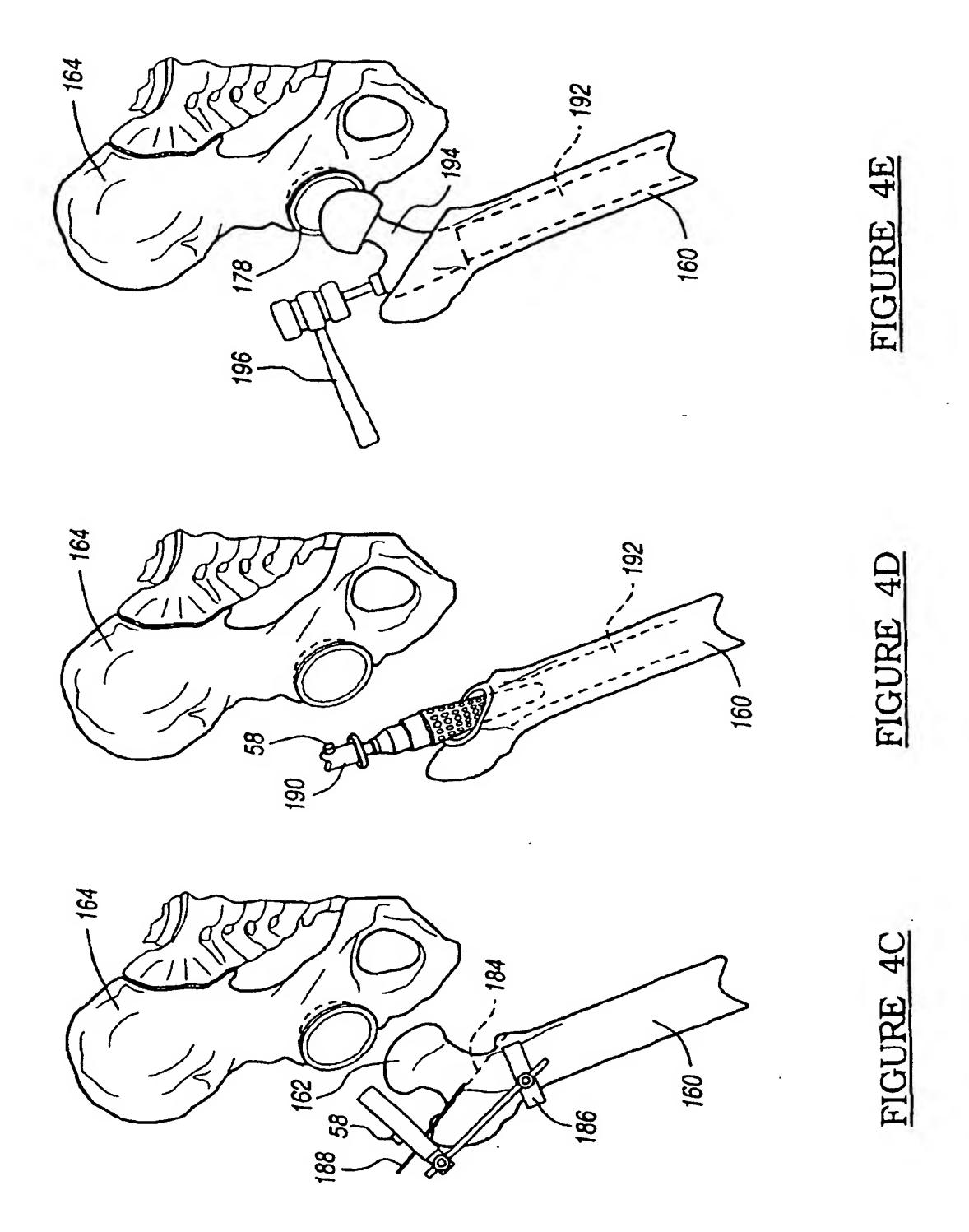


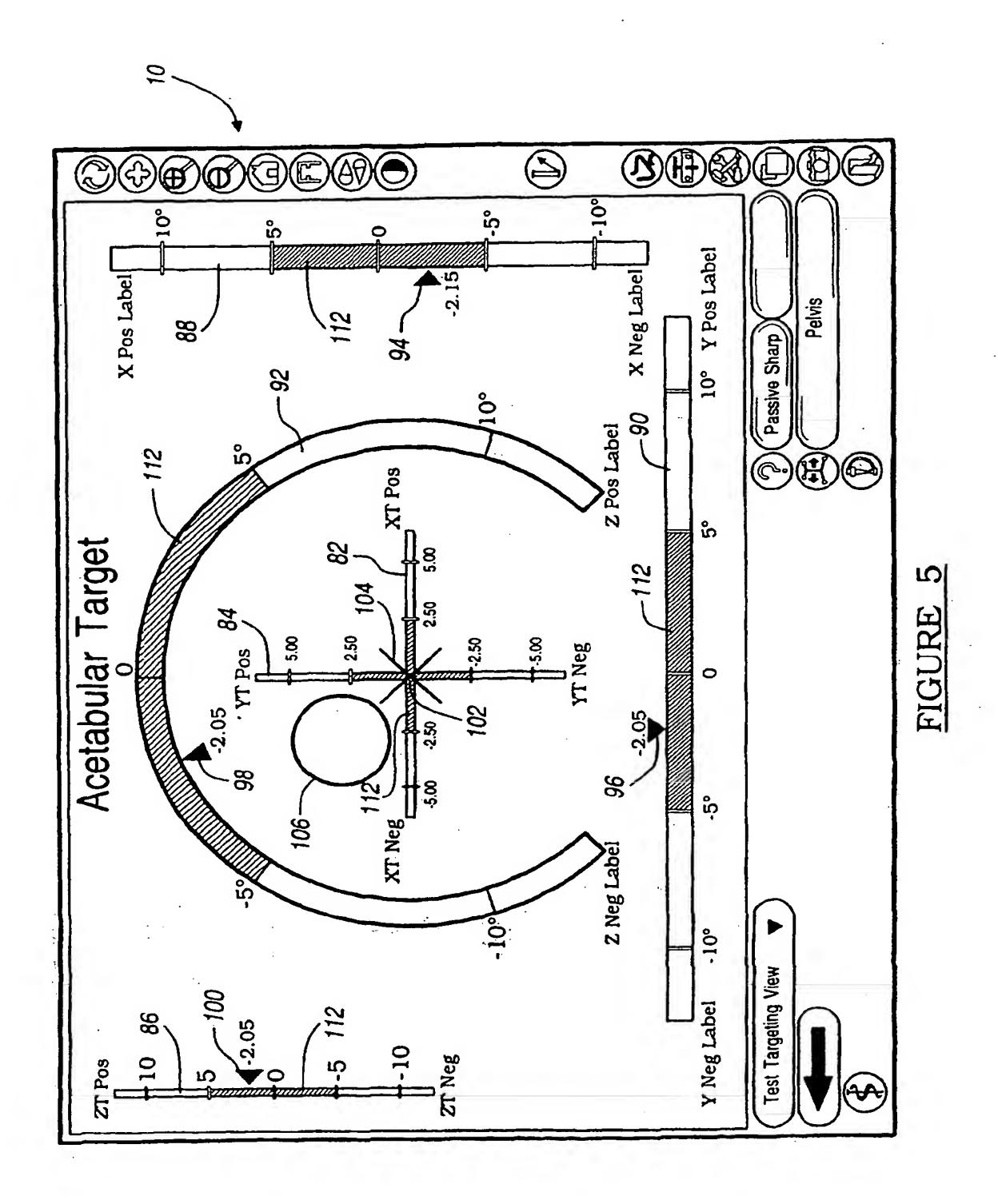

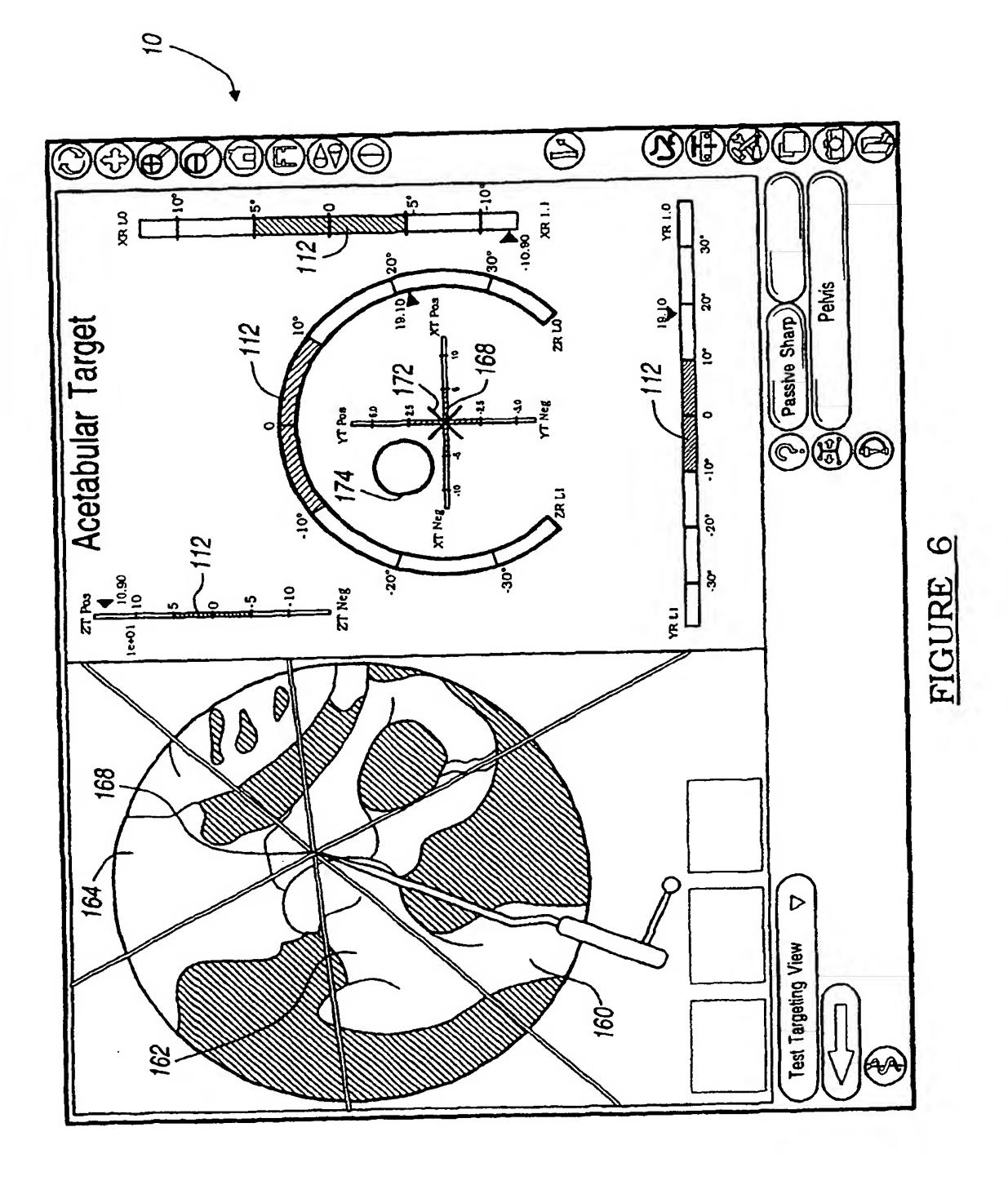

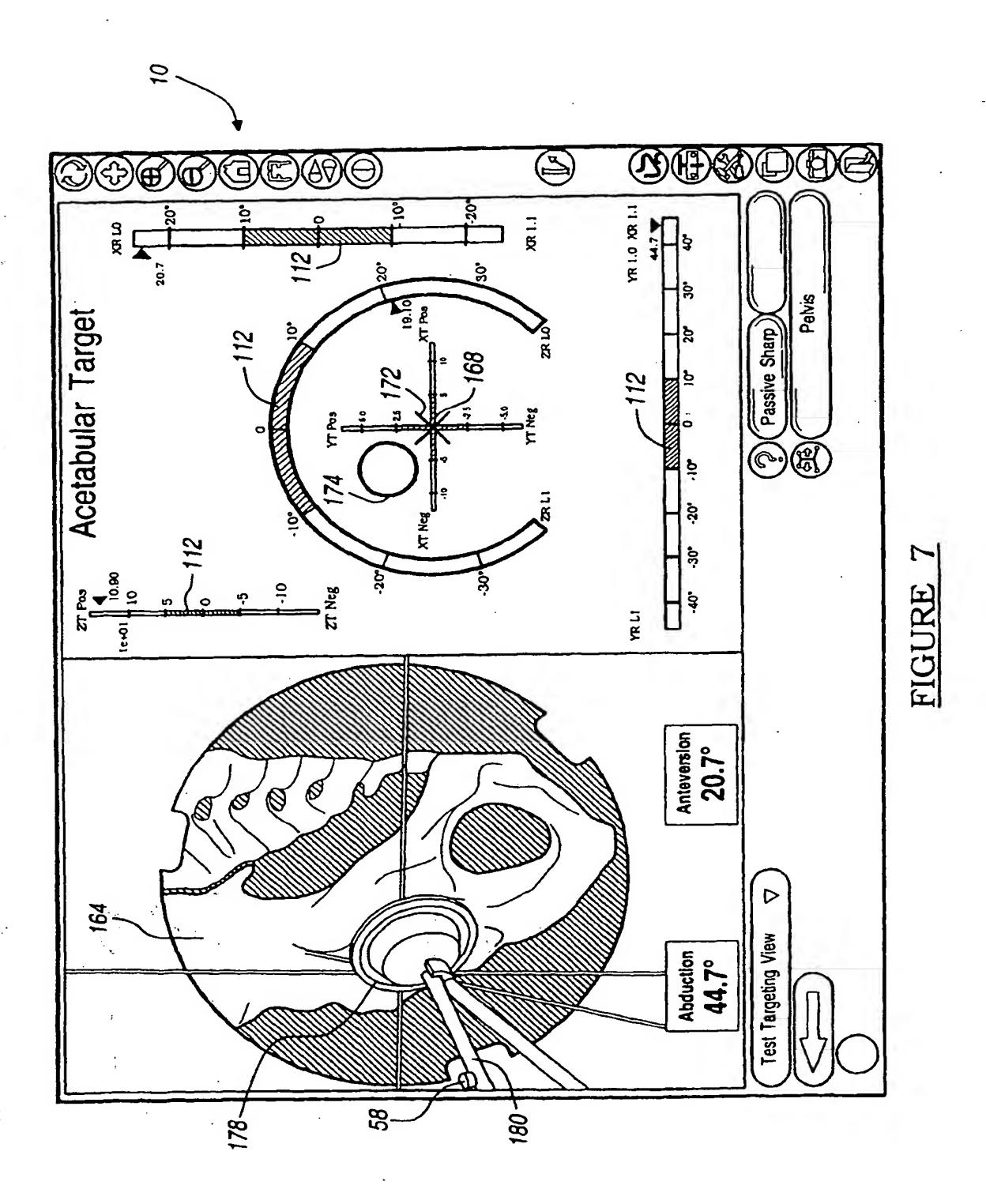

**3**5

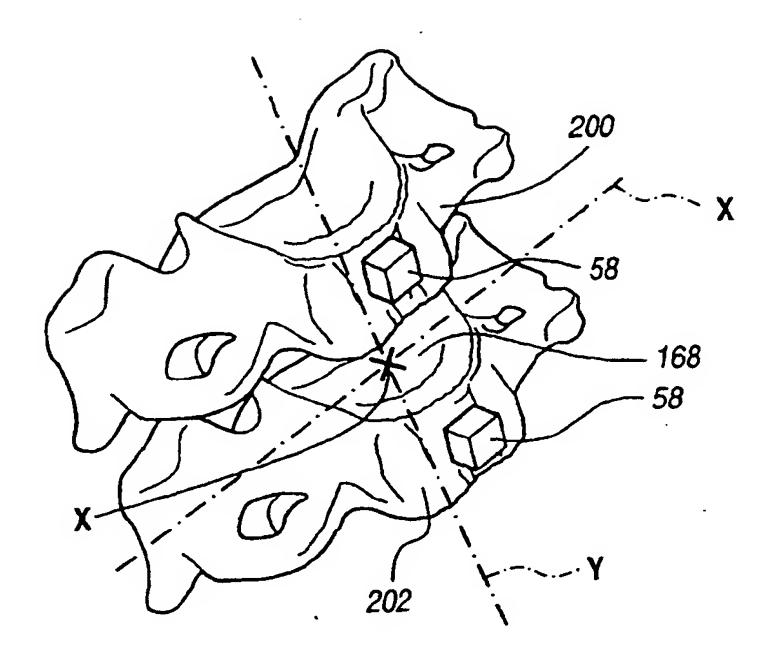



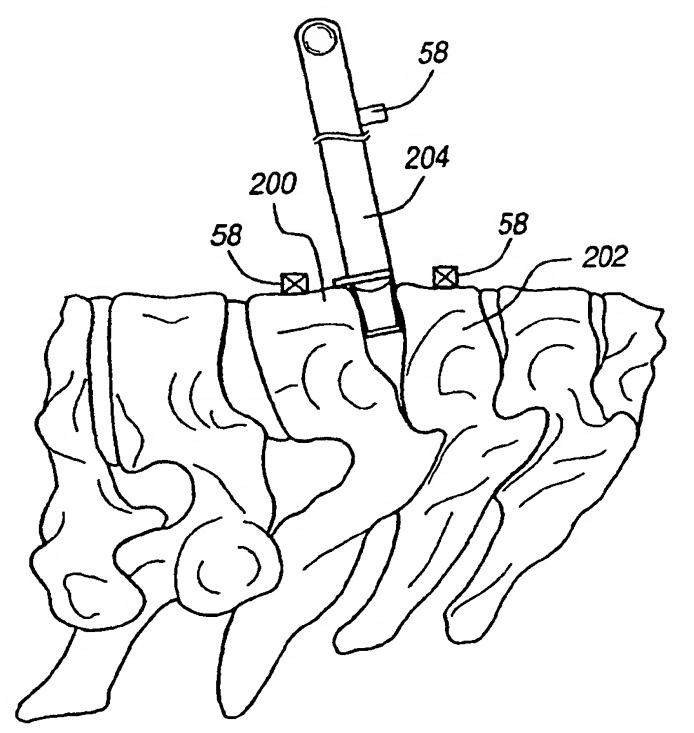
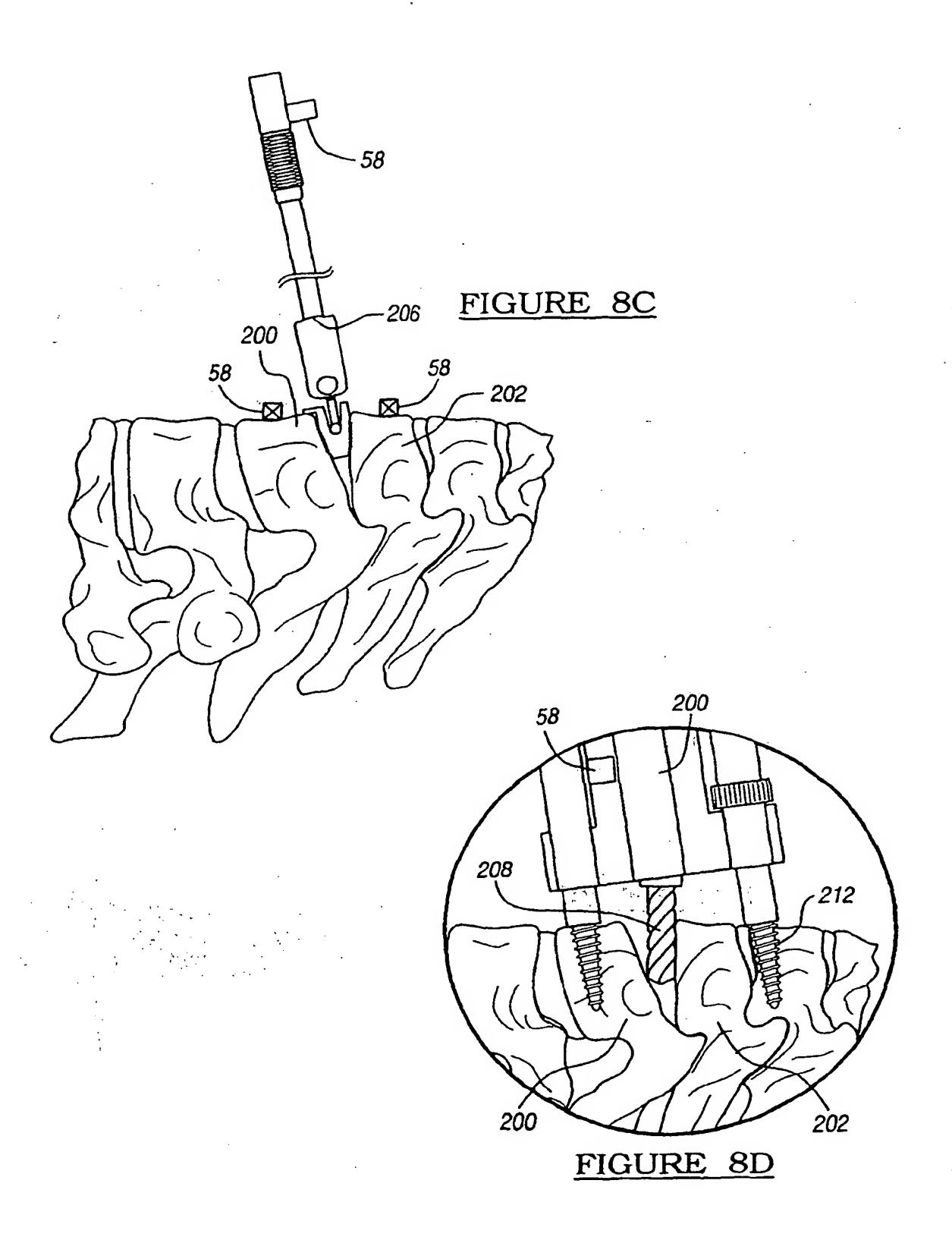


FIGURE 8B







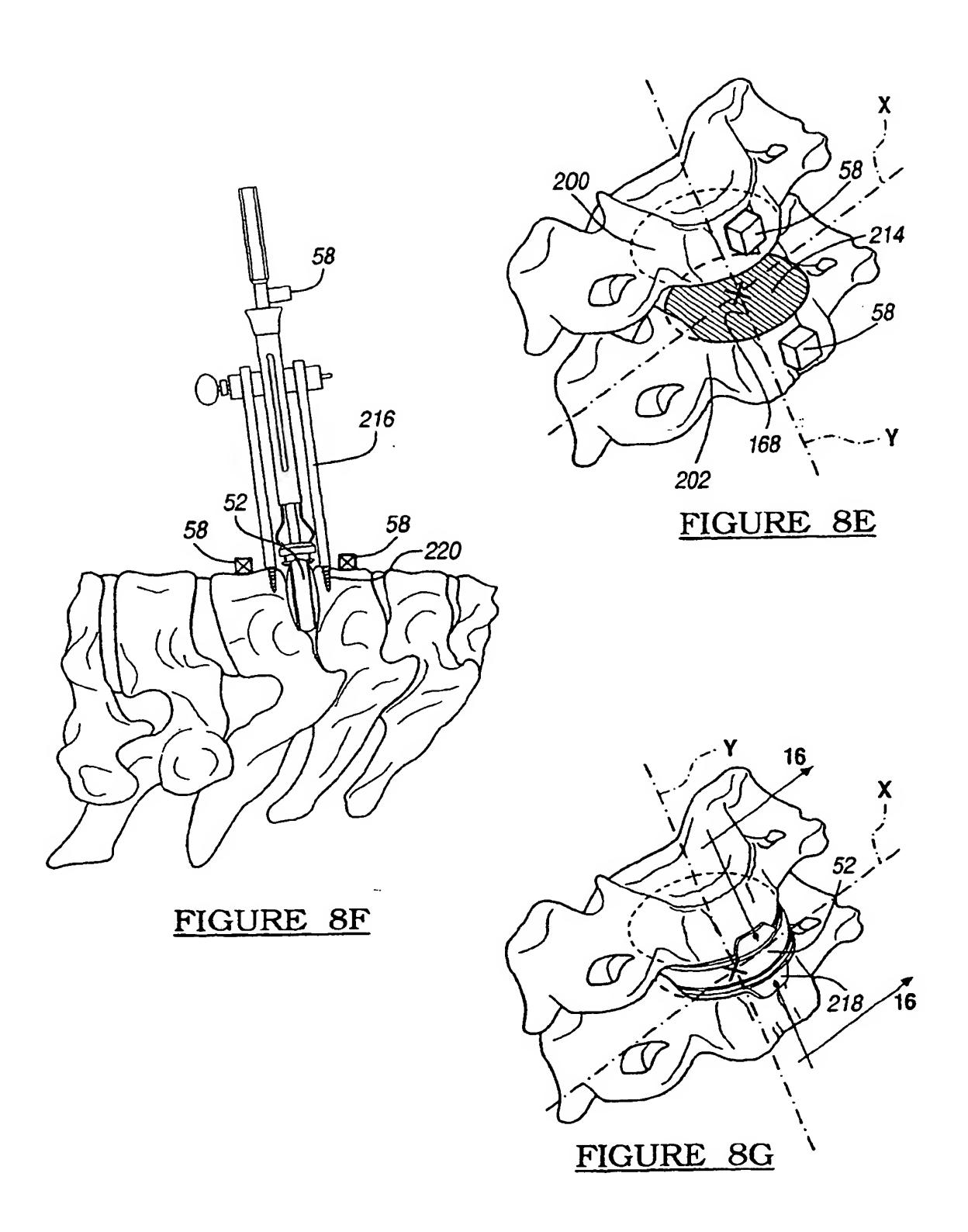



FIGURE 9

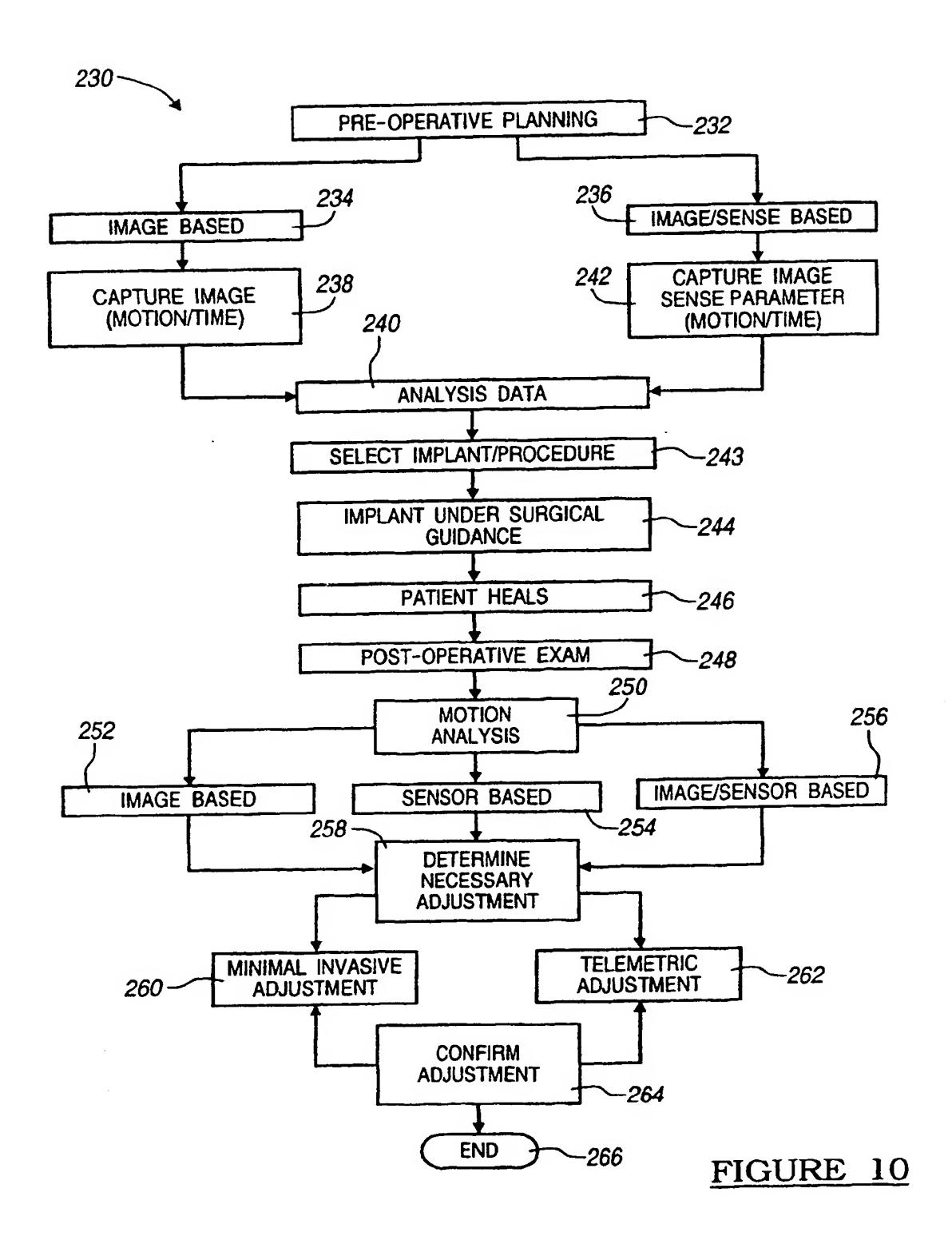


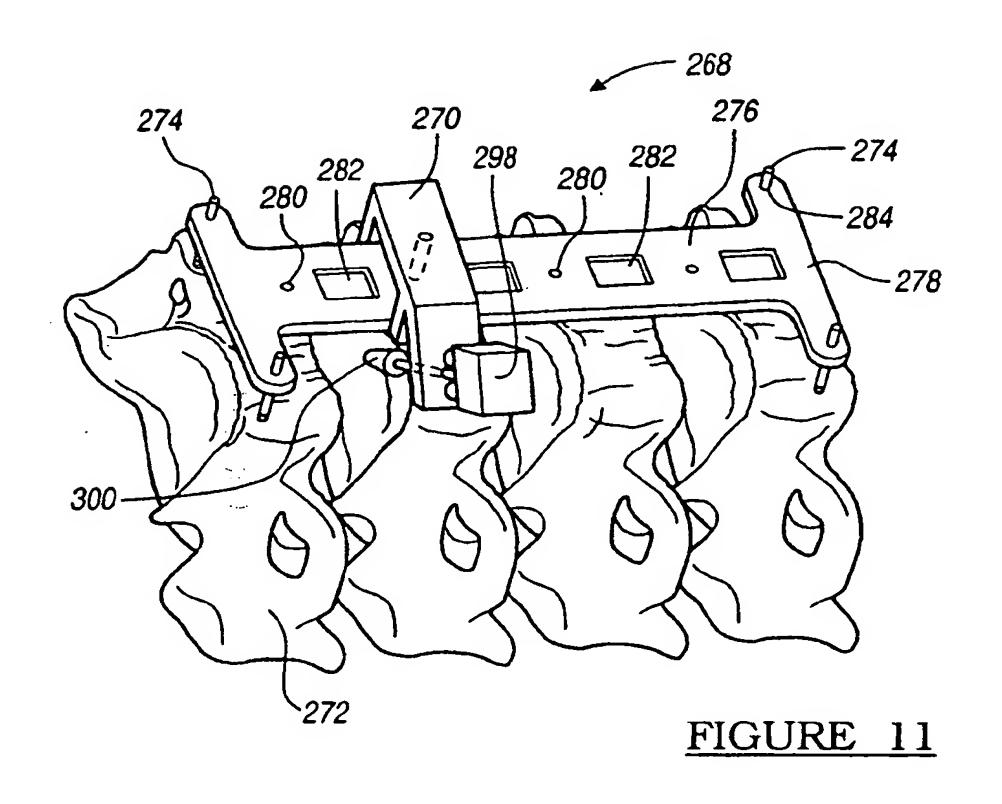


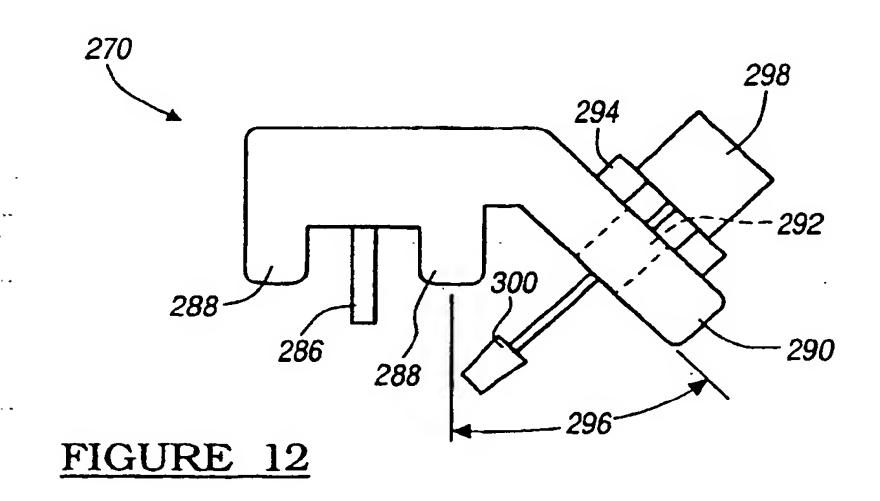


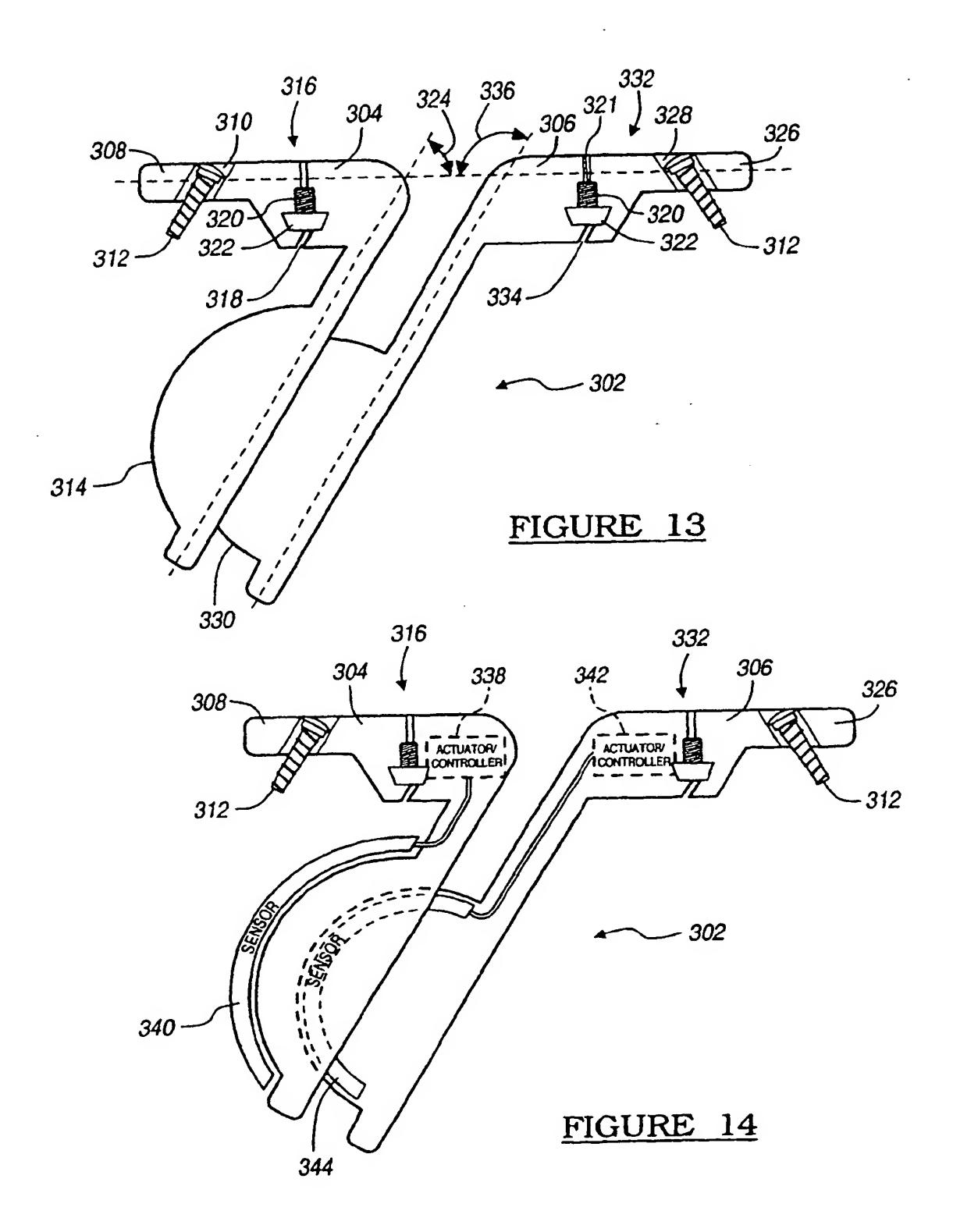


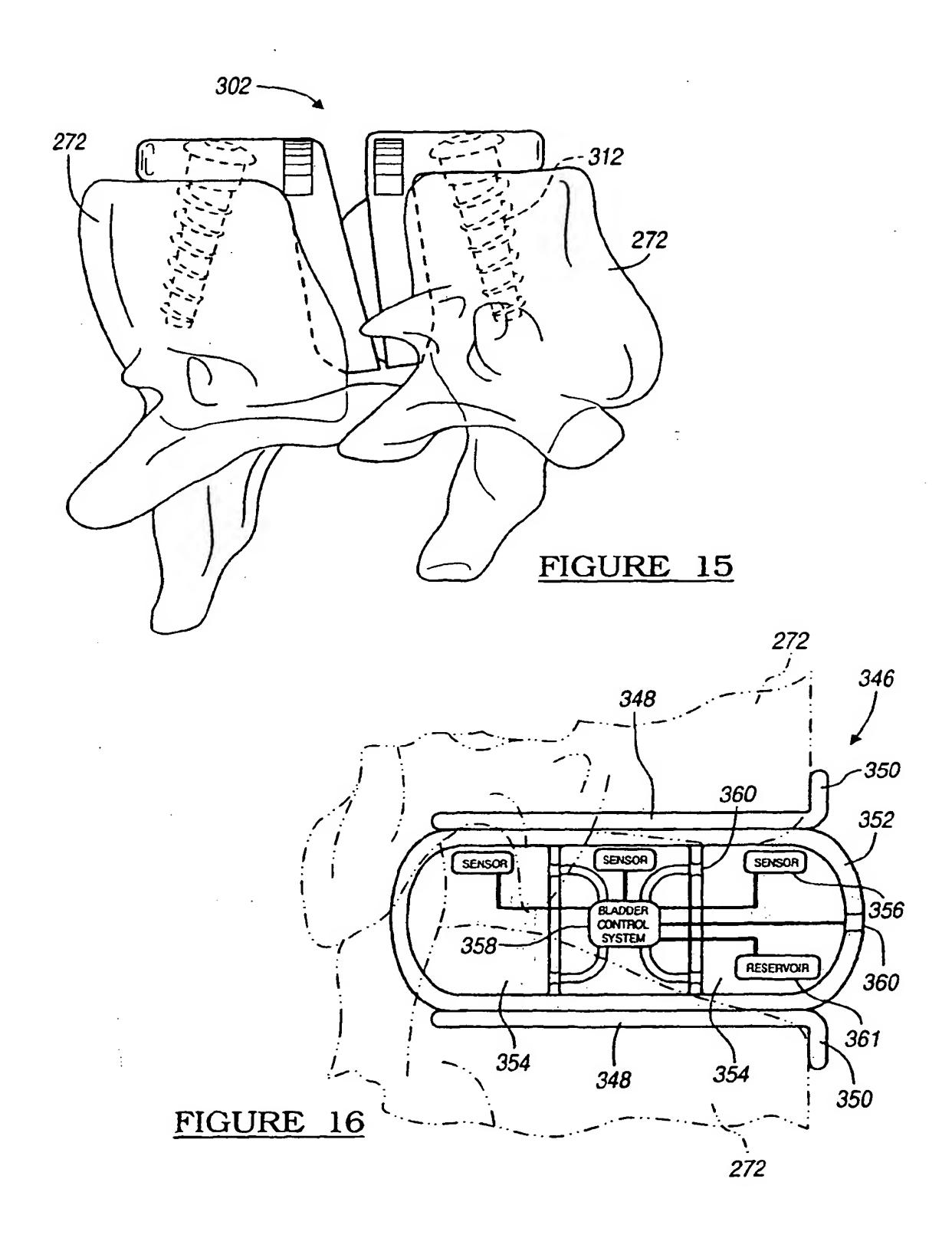


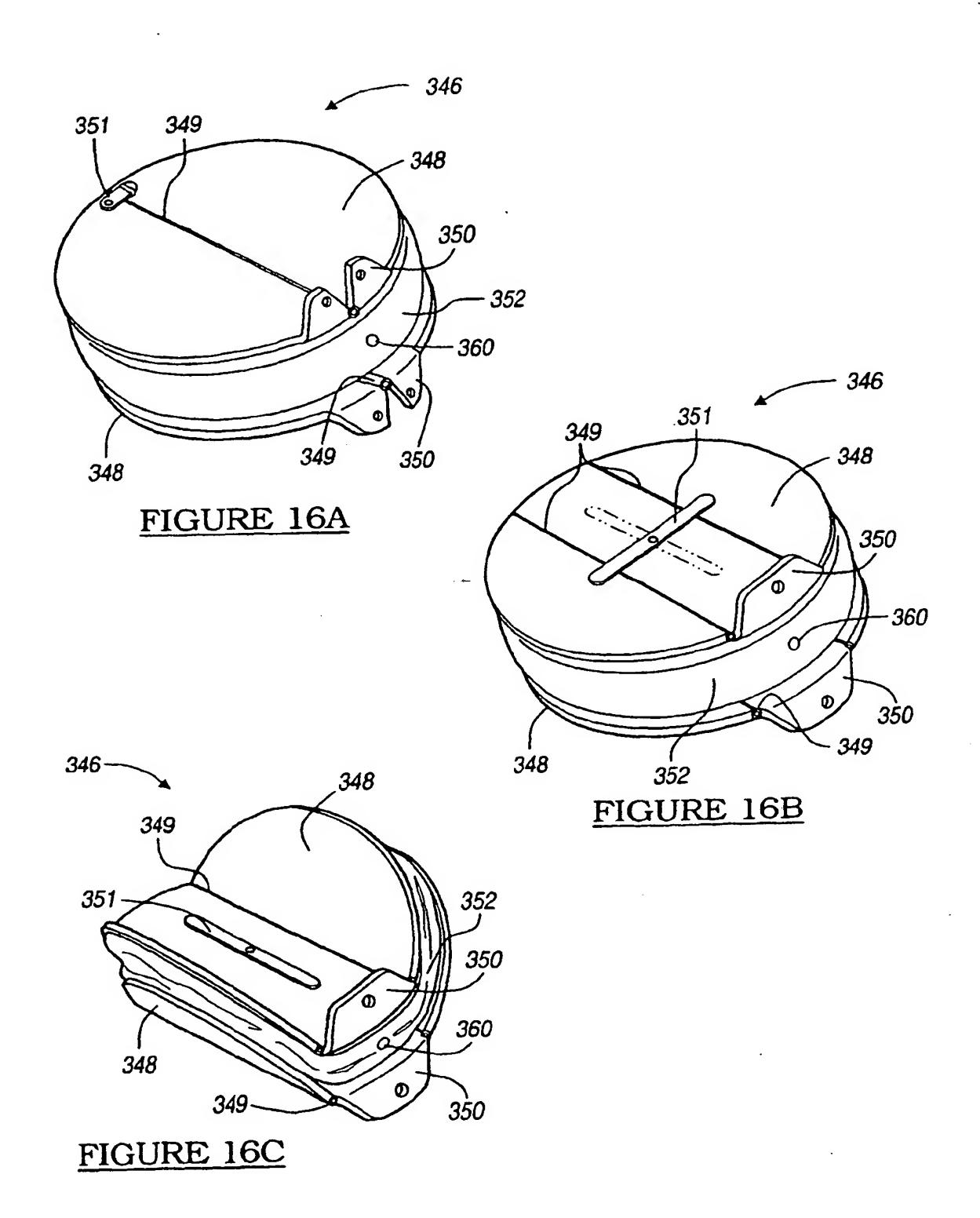

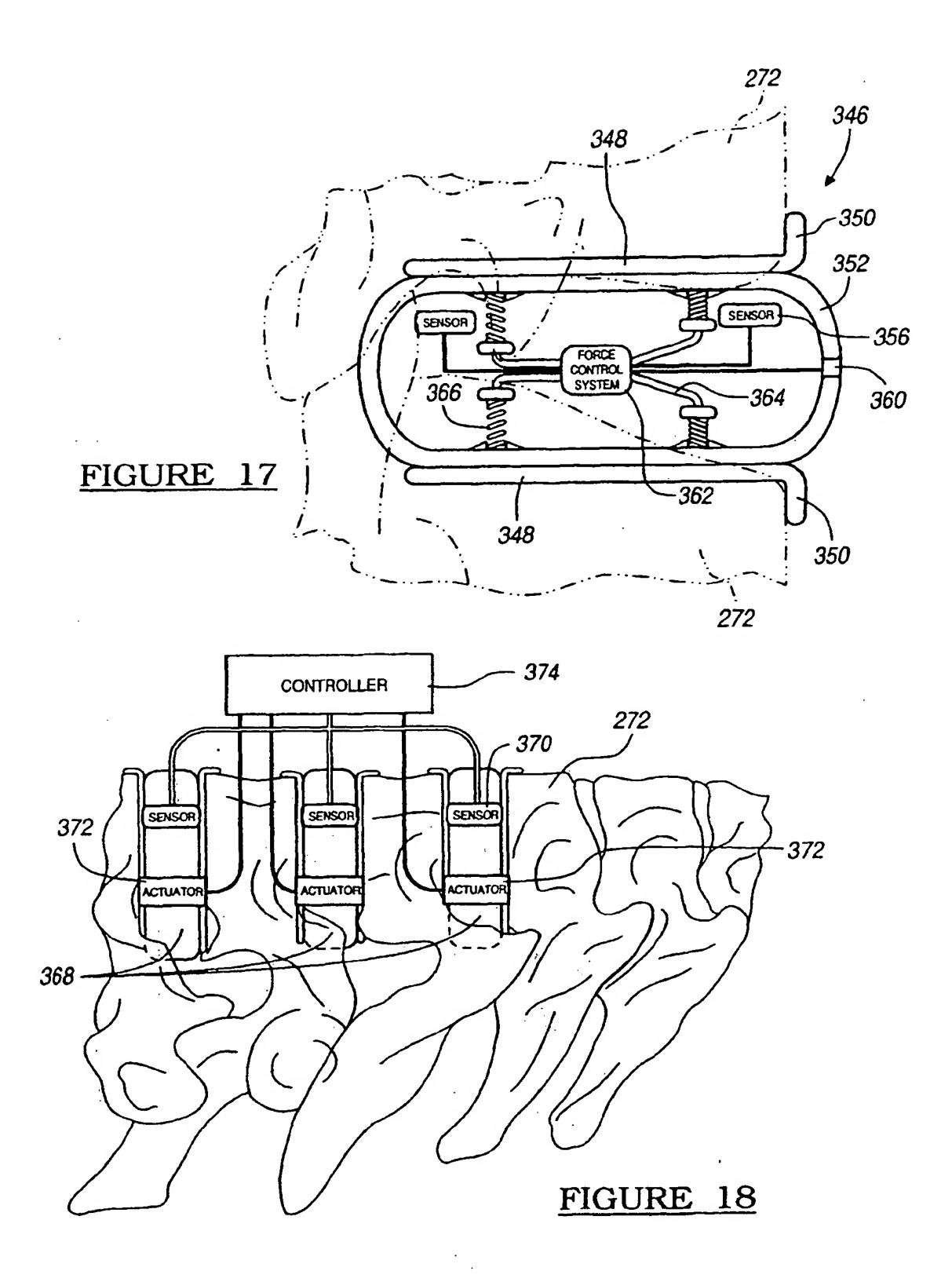

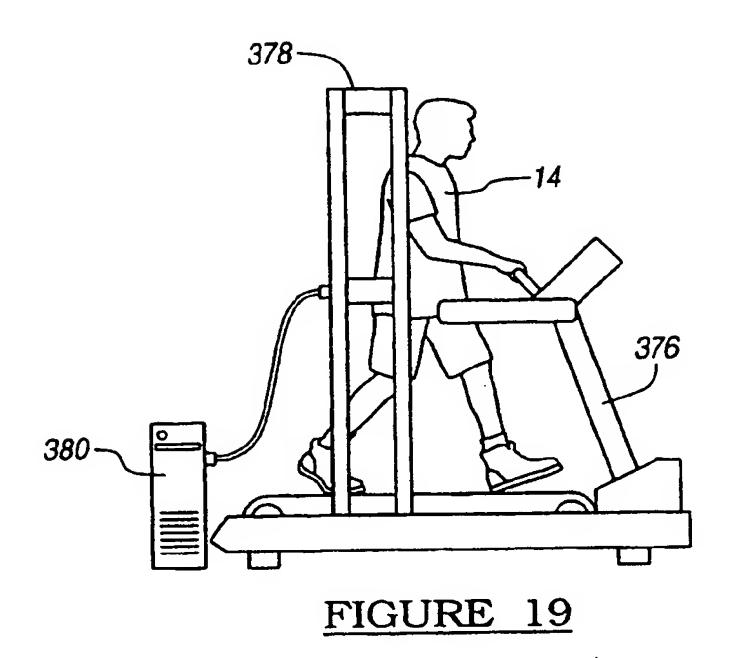

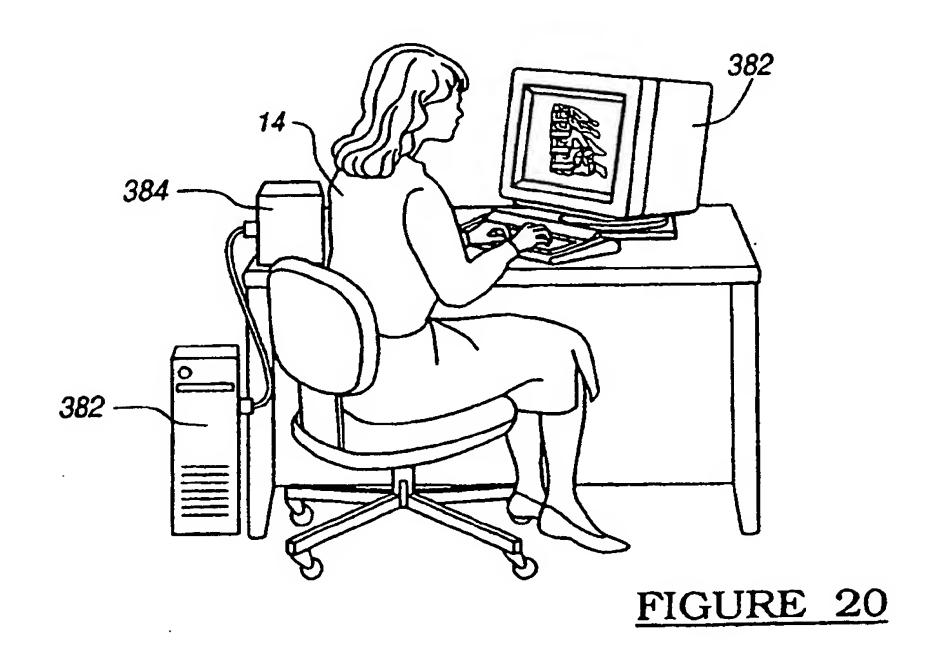

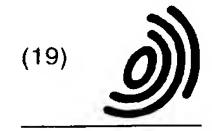

Europäisches Patentamt

European Patent Office

Office européen des brevets



(11) EP 1 442 715 A3

(12)

## **EUROPEAN PATENT APPLICATION**

- (88) Date of publication A3: 03.11.2004 Bulletin 2004/45
- (43) Date of publication A2: 04.08.2004 Bulletin 2004/32
- (21) Application number: 04001428.4
- (22) Date of filing: 23.01.2004

- (51) Int CI.7: **A61B 17/17**, A61B 19/00, A61F 2/32, A61F 2/44, A61F 2/46
- (84) Designated Contracting States:

  AT BE BG CH CY CZ DE DK EE ES FI FR GB GR

  HU IE IT LI LU MC NL PT RO SE SI SK TR

  Designated Extension States:


  AL LT LV MK
- (30) Priority: **30.01.2003 US 354562 25.04.2003 US 423515**

- (71) Applicant: SURGICAL NAVIGATION
  TECHNOLOGIES, INC.
  Louisville, Colorado 80027 (US)
- (72) Inventor: Hunter, Mark W. Broomfield Colorado 80020 (US)
- (74) Representative: Edlund, Fabian Awapatent AB,
  P.O. Box 11 394
  404 28 Göteborg (SE)

#### (54) Tunable spinal implant and apparatus for its post-operative tuning

(57) A tunable implant, system, and method enables a tunable implant (302) to be adjusted within a patient. The tunable implant includes a securing mechanism (308,310) to secure the implant in the patient, a actuation portion that enables the implant to move and an adjustment portion that permits adjustment of the implant after the implant has been positioned within the patient. The method of adjusting the tunable implant includes analyzing the operation of the implant, determining if

any adjustments are necessary and adjusting the implant to improve implant performance. The implant system includes both the tunable implant and a telemetric system (378,384) that is operable to telemetrically receive data from the tunable implant where the data is used to determine if adjustment of the tunable implant is necessary. The system also includes an instrument assembly that is used for performing spinal surgery where the instrument assembly includes a mounting platform (268) and a jig (270).

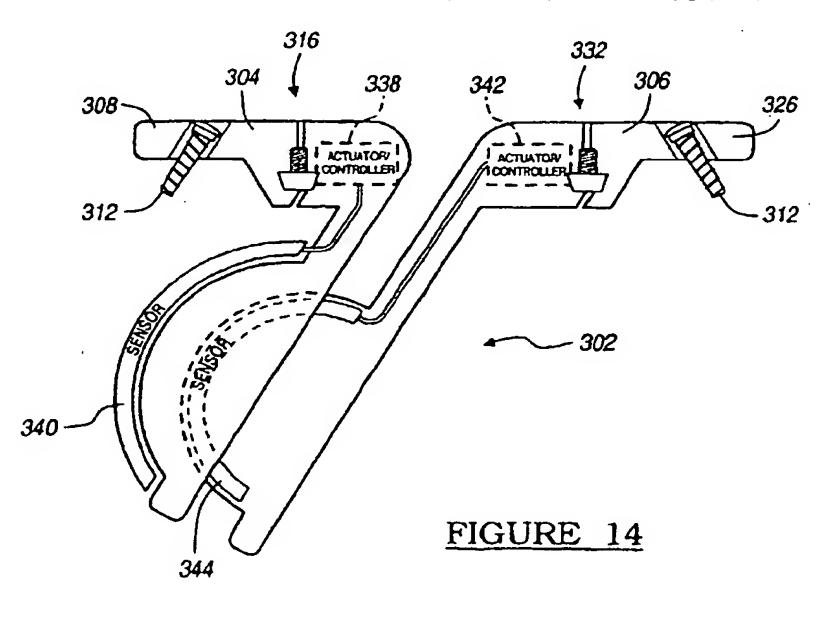

Printed by Jouve, 75001 PARIS (FR)



## PARTIAL EUROPEAN SEARCH REPORT

Application Number

which under Rule 45 of the European Patent Convention EP 04 00 1428 shall be considered, for the purposes of subsequent proceedings, as the European search report

	DOCUMENTS CONSIDE	RED TO BE RELEVANT			
ategory	Citation of document with income of relevant passage		Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.CI.7)	
(	EP 0 820 731 A (ZACC 28 January 1998 (199	98-01-28)	25-31, 33-39, 41,42	A61B17/17 A61B19/00 A61F2/32 A61F2/44	
	<pre>* column 7, line 18 * column 9, line 1 * column 15, line 29 * column 21, line 6</pre>	- column 6, line 20 * - line 58 * - line 31 * 5 - column 16, line 6 * - column 23, line 11 * 5 - column 27, line 27	-	A61F2/46	
(	US 6 190 414 B1 (YOU 20 February 2001 (20 * the whole document	901-02-20)	1,21,22		
<b>(</b>	GB 1 243 353 A (LUMI 18 August 1971 (1971 * page 3, line 63 -	1-08-18)	1-5		
<b>΄</b> Υ	l 8 June 2000 (2000-0	TTENSTEIN GMBH & CO KG) 6-08) - column 3, line 38 *	1,2,4-6, 23 24	A61B	
		-/		A61F	
INCO	MPLETE SEARCH		<u> </u>		
The Sear not comp be carrie	oh Division considers that the present	application, or one or more of its claims, does a meaningful search into the state of the art ca y, for these claims	/do annot		
Claims s	earched incompletely :				
Claims n	ot searched :				
Reason	for the limitation of the search:				
see	sheet C				
	Place of search	Date of completion of the search		Examiner	
	Munich	7 September 2004	Buc	chmann, G	
X:pa Y:pa	CATEGORY OF CITED DOCUMENTS rticularly relevant if taken alone rticularly relevant if combined with anot current of the same category	T : theory or principl E : earlier patent do after the filing da  her D : document cited i L : document cited f	cument, but publite in the application	ished on, or	
	chnological background n-written disclosure	same patent family, corresponding			



## INCOMPLETE SEARCH SHEET C

Application Number EP 04 00 1428

Claim(s) searched completely: 1-42 Claim(s) not searched: 53 Reason for the limitation of the search (non-patentable invention(s)): Article 52 (4) EPC - Method for treatment of the human or animal body by surgery



## PARTIAL EUROPEAN SEARCH REPORT

Application Number

EP 04 00 1428

	DOCUMENTS CONSIDERED TO BE RELEVANT	CLASSIFICATION OF THE APPLICATION (Int.CI.7)	
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	
X	DE 100 13 519 A (ADEVA MEDICAL GES FUER ENTWICK) 4 October 2001 (2001-10-04) * column 3, line 3 - column 4, line 47 *	1,7,8, 18-20,23	
X	US 4 584 994 A (BAMBERGER CHARLES ET AL) 29 April 1986 (1986-04-29) * the whole document *	1,7-11	
X Y	US 5 961 553 A (COTY ALAIN ET AL) 5 October 1999 (1999-10-05) * column 2, line 36 - column 6, line 48 *	33-39,41 40	
Y	US 5 415 660 A (BECHTOLD JOAN E ET AL) 16 May 1995 (1995-05-16) * column 7, line 23 - column 8, line 29 *	40	TECHNICAL FIELDS SEARCHED (Int.CI.7)
A	EP 1 057 461 A (WOLTER DIETMAR) 6 December 2000 (2000-12-06) * column 6, line 34 - line 46 *	33-42	
P,A	WO 03/079940 A (TRUSTEES OF THE UNIVERSITY OF) 2 October 2003 (2003-10-02) * paragraph [0033] - paragraph [0043] *	33-42	
A	WO 00/23015 A (SDGI HOLDINGS INC) 27 April 2000 (2000-04-27) * the whole document *	1-31	
A	US 5 865 846 A (BRYAN VINCENT ET AL) 2 February 1999 (1999-02-02) * the whole document *	1-31	
		i I	
		I	
İ	·	1	



**Application Number** 

EP 04 00 1428

CLAIMS INCURRING FEES
The present European patent application comprised at the time of filing more than ten claims.
Only part of the claims have been paid within the prescribed time limit. The present European search report has been drawn up for the first ten claims and for those claims for which claims fees have been paid, namely claim(s):
No claims fees have been paid within the prescribed time limit. The present European search report has been drawn up for the first ten claims.
LACK OF UNITY OF INVENTION
The Search Division considers that the present European patent application does not comply with the requirements of unity of invention and relates to several inventions or groups of inventions, namely:
see sheet B
All further search fees have been paid within the fixed time limit. The present European search report has been drawn up for all claims.
As all searchable claims could be searched without effort justifying an additional fee, the Search Division did not invite payment of any additional fee.
Only part of the further search fees have been paid within the fixed time limit. The present European search report has been drawn up for those parts of the European patent application which relate to the inventions in respect of which search fees have been paid, namely claims:
1-42
None of the further search fees have been paid within the fixed time limit. The present European search report has been drawn up for those parts of the European patent application which relate to the invention first mentioned in the claims, namely claims:



# LACK OF UNITY OF INVENTION SHEET B

**Application Number** 

EP 04 00 1428

The Search Division considers that the present European patent application does not comply with the requirements of unity of invention and relates to several inventions or groups of inventions, namely:

1. claims: 1-32

A tunable implant comprising a securing mechanism, an actuation portion and an adjustable portion.

2. claims: 33-42

A tunable implant system comprising an implant with an adjustment portion, and a telemetric system for adjustment of the implant.

3. claims: 43-52

An instrument assembly for use in spinal surgery, comprising a mounting platform to be positioned adjacent to a plurality of vertebrae, and a jig to be attached to the mounting platform and to support an instrument.

## ANNEX TO THE EUROPEAN SEARCH REPORT ON EUROPEAN PATENT APPLICATION NO.

EP 04 00 1428

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report. The members are as contained in the European Patent Office EDP file on The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

07-09-2004

	Patent document ted in search report		Publication date		Patent family member(s)	•	Publicatio date
EP	0820731	A	28-01-1998	FR FR AU CA DE EP US	2751201 2751202 732244 2873497 2213063 69722320 0820731 2002151978	A1 B2 A A1 D1 A2	23-01-1 23-01-1 12-04-2 29-01-1 22-01-1 03-07-2 28-01-1 17-10-2
US	6190414	B1	20-02-2001	NONE			
GB	1243353	Α	18-08-1971	DE US	1807227 3426364		27-05-1 11-02-1
DE	19856013	Α	08-06-2000	DE WO EP JP	19856013 0033752 1135076 2002531214	A1 A1	08-06-2 15-06-2 26-09-2 24-09-2
DE	10013519	Α	04-10-2001	DE WO EP	10013519 0170131 1265551	A1	04-10-20 27-09-20 18-12-20
US	4584994	Α	29-04-1986	NONE			~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~ ~
US	5961553	A	05-10-1999	FR AT DE DE EP WO	2730406 172373 69600837 69600837 0809473 9625117	T D1 T2 A1	14-08-19 15-11-19 26-11-19 27-05-19 03-12-19 22-08-19
US	5415660	Α	16-05-1995	NONE			
EP .	1057461	A	06-12-2000	DE US CA EP JP	19924676 6706071 2307993 1057461 2001046411	B1 A1 A1	30-11-20 16-03-20 29-11-20 06-12-20 20-02-20
WO	03079940	A	02-10-2003	WO US	03079940 2004019382	•	02-10-20 29-01-20
WO	0023015	Α	27-04-2000	WO AU EP	0023015 1116799 1124508	Α	27-04-20 08-05-20 22-08-20

## ANNEX TO THE EUROPEAN SEARCH REPORT ON EUROPEAN PATENT APPLICATION NO.

EP 04 00 1428

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report. The members are as contained in the European Patent Office EDP file on The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

07-09-2004

			Publication date		Patent family member(s)	Publication date
WO	0023015	A		JP	2002508679 T	19-03-2002
US	5865846	A	02-02-1999	CA US EP US AT DE EP ES US	2202453 A1 5674296 A 0820740 A1 6001130 A 253339 T 69725932 D1 69725932 T2 1166725 A2 2210458 T3 6156067 A	11-10-1998 07-10-1998 28-01-1998 14-12-1999 15-11-2003 11-12-2003 02-09-2004 02-01-2003 01-07-2004
<b></b> .	<b></b>		, <b>, , , , , , , , , , , , , , , , , , </b>			
					•	
		٠				
	WO	Palent document cited in search report  WO 0023015  US 5865846	wo 0023015 A	cited in search report date  W0 0023015 A	Cited in search report date  WO 0023015 A JP  US 5865846 A 02-02-1999 CA	cited in search report         date         member(s)           W0 0023015         A         JP 2002508679 T           US 5865846         A 02-02-1999         CA 2202453 A1           US 5674296 A         EP 0820740 A1           US 6001130 A         AT 253339 T           DE 69725932 D1         DE 69725932 T2           EP 1166725 A2         ES 2210458 T3